

पारम्परिक औषधियों में उपयोग हेतु अंजीर  
दशती (*Ficus hispida* L. f.) सूखा  
फल — विशिष्टि

**Anjīr Dashtī (*Ficus hispida* L. f.)  
Dried Fruit for Use in Traditional  
Medicine — Specification**

ICS 11.120.10

© BIS 2024



भारतीय मानक ब्यूरो  
BUREAU OF INDIAN STANDARDS  
मानक भवन, 9 बहादुर शाह ज़फर मार्ग, नई दिल्ली - 110002  
MANAK BHAVAN, 9 BAHADUR SHAH ZAFAR MARG  
NEW DELHI - 110002  
[www.bis.gov.in](http://www.bis.gov.in) [www.standardsbis.in](http://www.standardsbis.in)

July 2024

Price Group 9

## FOREWORD

This Indian Standard was adopted by Bureau of Indian Standards, after the draft finalized by *Unani* Sectional Committee had been approved by the Ayush Division Council.

*Anjīr dashī* a moderate sized tree or shrub, distributed throughout the outer *Himalayan* range from *Chenab* eastwards to Bengal, Central and South India and Andaman Islands. *Anjīr dashī* consists of dried fruits of *Ficus hispida* L.f. (Family Moraceae). The dried fruits of *Anjīr Dashī* used as an ingredient of different formulations (*Mā' al-Lahm Khās*, *Safoof Bars*, *La 'ūq-i Dīq al-Nafas Balghamī*) in Unani medicine .

*Anjīr Dashī* is known by different names in other languages such as *Tīn* (Arabic); *Khoskadumar*, *Tanvardi*, *Teenbarree* (Assamese); *Kakdumur*, *Kathdumur*, *Kakadumbar* (Bengali); Wild Fig, Devil Fig (English); *Tedumbaro*, *Dhedadambaro*, *Dhedhumbro* (Gujarati); *Konea-dumbar*, *Kathumar* (Kannad); *Kadaatti*, *Arjeeru Hamu*, *Anjeeru*, *Onagida*, *Hanna* (Hindi); *Peyatti*, *Kattatti*, *Erumanakku*, *Parakasimi* (Malayalam); *Rambal*, *Kalodumbar*, *Bhuiumbar* (Marathi); *Dimiri*, *Ani Dambura* (Oriya); *Anjīr Dashī* (Persian); *Rumbal* (Punjabi); *Kakodumbur*, *Malayu*, *Malpu* (Sanskrit); *Peyatti* (Tamil); *Brahmamedi*, *Kakimedi* (Telugu); *Kathguular* (Urdu).

This standard is one of the series of standards being brought out on ingredients used in formulations of Unani medicine for the advantage of researchers, academicians, students, clinical practitioners and drug manufacturers, etc.

In the formulation of this standard, significant assistance has been derived from the *Unani Pharmacopoeia of India*, Part 1, Vol VII, 2022 published by the Pharmacopoeia Commission for Indian Medicine and Homoeopathy, Ministry of Ayush, Government of India. Inputs have also been derived from the information available in the public domain in print and electronic media including authoritative books.

In the formulation of this standard due consideration has been given to the provisions of *The Drugs and Cosmetics Act*, 1940 and Rules framed thereunder. However, this standard is subject to the restrictions imposed under these Rules and Regulations, wherever applicable.

The composition of the Committee responsible for the formulation of this standard is given in [Annex D](#).

For the purpose of deciding whether a particular requirement of this standard is complied with, the final value, observed or calculated, expressing the result of a test or analysis, shall be rounded off in accordance with IS 2 : 2022 'Rules for rounding off numerical values (*second revision*)'. The number of significant places retained in the rounded off value shall be the same as that of the specified value in this standard.

*Indian Standard*

## ANJĪR DASHTĪ (*Ficus hispida* L. f.) DRIED FRUIT FOR USE IN TRADITIONAL MEDICINE — SPECIFICATION

### 1 SCOPE

This standard prescribes the requirements and methods of test for *Anjīr dashtī* which consists of dried fruits of (*Ficus hispida* L.f.), Family of Moraceae.

### 2 REFERENCES

The standards listed in [Annex A](#) contain provisions which, through reference in this text, constitute provisions of this standard. At the time of publication, the editions indicated were valid. All standards are subject to revision and parties to agreements based on this standard are encouraged to investigate the possibility of applying the most recent edition of these standards.

### 3 REQUIREMENTS

#### 3.1 Description

##### 3.1.1 Macroscopic Examination of *Anjīr dashtī* Dried Fruit

*Anjīr dashtī* is dried syconus ovoid fruit with a central circular hole and short stalk. The fruit is 1 cm to 2 cm in diameter, greyish brown in colour and wrinkled. Seeds are less than 1 mm in diameter, yellowish-brown in colour. Odour and taste are not characteristic.

##### 3.1.2 Microscopic Examination of *Anjīr dashtī* Dried Fruit

Fruit shows a single layered epidermis, covered with thick cuticle and few unicellular trichomes. Epidermis are 4 to 6 layers of hexagonal to polygonal collenchymatous cells. A few cells are containing rosette crystals of calcium oxalate. Mesocarp is composed of large, oval to polygonal, thick-walled parenchymatous cells. A few vascular vessels are showing spiral thickening.

*Anjīr dashtī powder* is brownish yellow in colour. Odour is not characteristic. Powder taste is astringent. Epidermal cells, unicellular trichomes and anomocytic stomata are present. Trichomes are

simple and glandular upto 350  $\mu$ . Simple trichomes are unicellular short, straight or bent and warty and long thick walled along with bulbous base. Rosette of calcium oxalate crystals measuring upto 40  $\mu$ , are found scattered or embedded in the collenchymatous cells of the receptacle. Spiral vessels are upto 25  $\mu$  along with fragments of septate laticifers vessels. Thick walled pitted sclereids are upto 100 in length. Epidermal cells of the seed are coated with endosperm cells and parenchyma cells (see [Fig. 1](#)).

##### 3.1.3 Macroscopic and Microscopic Examination of *Anjīr Dashtī* Powder

#### 3.2 General

**3.2.1** *Anjīr dashtī* dried fruits shall be free from extraneous/artificial flavours.

**3.2.2** *Anjīr dashtī* dried fruits shall comply with physical, chemical and microbiological requirements given in [Table 1](#).

**3.2.3** If any pesticides other than those for which minimum requirements are given in [Table 1](#) are applied to the herb before or after harvesting, those should also be tested. Their limit shall be calculated using the following formula:

$$\frac{ADI \times M}{MDD \times 100}$$

where

ADI = acceptable daily intake of pesticide as published by [Food and Agriculture Organization-World Health Organization (FAO-WHO)], in milligrams per kilogram of body mass;

*M* = body mass, in, kilograms (60 kg); and

*MDD* = maximum daily dose of the drug, in kilograms.

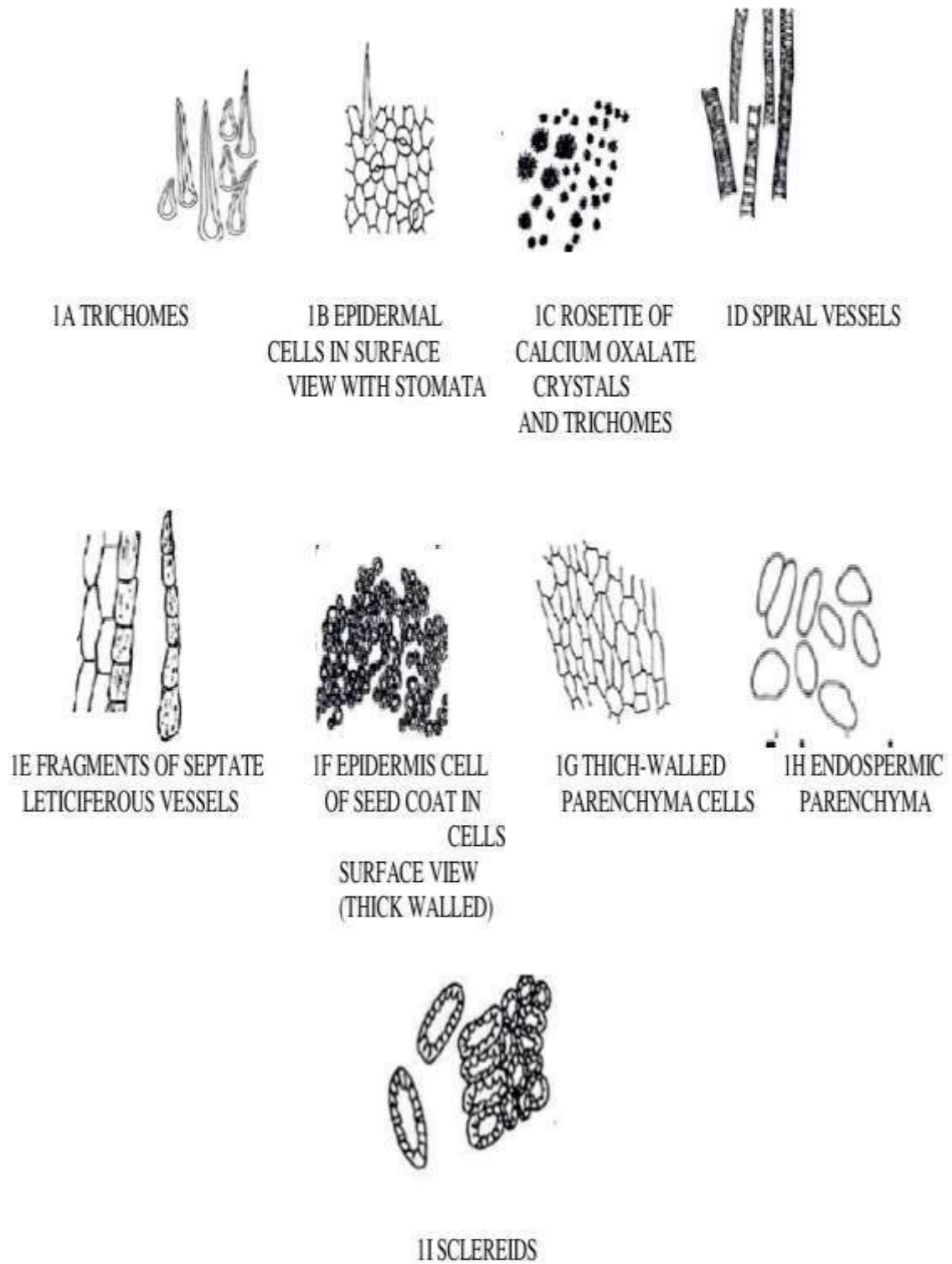


FIG. 1 ANJĪR DASHTĪ — POWDER MICROSCOPY

## 4 PACKING, STORAGE AND MARKING

### 4.1 Packing

*Anjīr dashī* dried fruits shall be packed in clean, sound and dry container made up of metal, glass, food grade polymers, wood or jute bags. The wooden boxes or jute bags shall be suitably lined with moisture proof lining which does not impart any foreign smell to the product. The packing material shall be free from any fungal or insect infestation and should not impart any foreign smell. Each container shall be securely closed and sealed.

### 4.2 Storage

*Anjīr dashī* dried fruits shall be stored under conditions that prevent contamination and, as far as possible, deterioration. Storage area shall be clean, well ventilated, protected from light, moisture, insects and rodents.

### 4.3 Marking

The following particulars shall be legibly and indelibly marked or labelled on each pack:

- Name of the material including part of the plant, botanical name, and trade name or brand name, if any;
- Name and address of the producer or packer;
- Batch number;
- Net quantity;
- State and country of production;
- Date of packing (MM/YYYY);
- Instructions for storage; and

- Any other information requested by the buyer, such as the date of harvesting (MM/YYYY) (if known).

The above information, or part of it, may instead appear in the documentation after agreement between the buyer and the seller.

### 4.4 BIS Certification Marking

The product(s) conforming to the requirements of this standard may be certified as per the conformity assessment schemes under the provisions of the *Bureau of Indian Standards Act, 2016* and the Rules and Regulations framed thereunder, and the product(s) may be marked with the Standard Mark.

## 5 SAMPLING

**5.1** Representative samples of the material shall be drawn and tested for conformity to this specification as prescribed in IS 13145.

**5.2** The samples of *Anjīr dashī* dried fruits shall be tested for ascertaining conformity of the material to the requirements in accordance with the relevant clauses given in col (4) of [Table 1](#).

## 6 QUALITY OF REAGENTS

**6.1** Reagents including pure chemicals used shall be of analytical grade.

**6.2** Reagent grade water for laboratory use shall be as per IS 1070.

NOTE — 'Pure chemicals' shall mean chemicals that do not contain impurities which effect the results of analysis.

**Table 1 Requirements for *Anjīr Dashī* (*Ficus hispida* L. f.) Dried Fruits**

(Clauses [3.2.2](#), [3.2.3](#) and [5.2](#))

SI No.	Characteristic	Requirement	Method of Test, Ref to
(1)	(2)	(3)	(4)
i)	Foreign matter, percent by mass, <i>Max</i>	1.0	<b>6.2</b> of IS 4333 (Part 1)
ii)	Total ash, percent by mass, <i>Max</i>	7.0	<b>6</b> of IS 1797
iii)	Acid insoluble ash, percent by mass, <i>Max</i>	1.5	<b>8</b> of IS 1797
iv)	Alcohol soluble extractive, percent by mass, <i>Min</i>	3.0	<b>10</b> of IS 1797

Table 1 (Continued)

SI No.	Characteristic	Requirement	Method of Test, Ref to
(1)	(2)	(3)	(4)
v)	Water soluble extractive, percent by mass, <i>Min</i>	12.0	11 of IS 1797
vi)	Volatile oil, percent by mass, <i>Min</i>	0.3	15 of IS 1797
vii)	Thin-layer chromatogram identification	-	<a href="#">Annex B</a>
viii)	Aflatoxin (B <sub>1</sub> ), ng/g, <i>Max</i>	2.0	IS 16287
ix)	Aflatoxin (B <sub>1</sub> + B <sub>2</sub> + G <sub>1</sub> + G <sub>2</sub> ), ng/g, <i>Max</i>	5.0	IS 16287
x)	Lead, mg/kg, <i>Max</i>	10.0	IS 16913
xi)	Mercury, mg/kg, <i>Max</i>	1.0	IS 16913
xii)	Cadmium, mg/kg, <i>Max</i>	0.3	IS 16913
xiii)	Arsenic, mg/kg, <i>Max</i>	3.0	IS 16913
xiv)	<i>Staphylococcus aureus</i> , per g	Absent for extract and powder Absent for plant material	<a href="#">Annex C</a>
xv)	<i>Pseudomonas aeruginosa</i> , per g	Absent for extract and powder Absent for plant material	<a href="#">Annex C</a>
xvi)	<i>Salmonella</i> spp., per g	Absent	<a href="#">Annex C</a>
xvii)	<i>Escherichia coli</i> , per g	Absent for extract and powder 10 for plant material	<a href="#">Annex C</a>
xviii)	Total microbial plate count, per g, <i>Max</i>	10 <sup>5</sup> for extract and powder 10 <sup>7</sup> for plant material	<a href="#">Annex C</a>
xix)	Total yeast and mould, per g, <i>Max</i>	10 <sup>3</sup> for extract and powder 10 <sup>5</sup> for plant material	<a href="#">Annex C</a>
xx)	Pesticide residues		
a)	Alachlor, mg/kg, <i>Max</i>	0.02	IS 17924
b)	Aldrin and dieldrin (sum of ), mg/kg, <i>Max</i>	0.05	
c)	Azinphos-methyl, mg/kg, <i>Max</i>	1.0	
d)	Bromopropylate, mg/kg, <i>Max</i>	3.0	
e)	Chlordane (sum of cis-, trans- and oxythlordane), mg/kg, <i>Max</i>	0.05	
f)	Chlorfenvinphos, mg/kg, <i>Max</i>	0.5	
g)	Chlorpyrifos, mg/kg, <i>Max</i>	0.2	

Table 1 (Concluded)

SI No.	Characteristic	Requirement	Method of Test, Ref to
(1)	(2)	(3)	(4)
h)	Chlorpyrifos-methyl, mg/kg, <i>Max</i>	0.1	IS 17924
j)	Cypermethrin (and isomers), mg/kg, <i>Max</i>	1.0	
k)	DDT (sum of p,p'-DDT, o,p'-DDT, p,p-DDE and p,p'-TDE), mg/kg, <i>Max</i>	1.0	
m)	Deltamethrin, mg/kg, <i>Max</i>	0.5	
n)	Diazinon, mg/kg, <i>Max</i>	0.5	
p)	Dichlorvos, mg/kg, <i>Max</i>	1.0	
q)	Dithiocarbamates (as CS <sub>2</sub> ), mg/kg, <i>Max</i>	2.0	
r)	Endosulfan (sum of isomers and endosulfan sulphate), mg/kg, <i>Max</i>	3.0	
s)	Endrin, mg/kg, <i>Max</i>	0.05	
t)	Ethion, mg/kg, <i>Max</i>	2.0	
u)	Fenitrothion, mg/kg, <i>Max</i>	0.5	
v)	Fenvalerate, mg/kg, <i>Max</i>	1.5	
w)	Fonofos, <i>Max</i>	0.05	
y)	Heptachlor (sum of heptachlor and heptachlor epoxide), mg/kg, <i>Max</i>	0.05	
z)	Hexachlorobenzene, mg/kg, <i>Max</i>	0.1	
aa)	Hexachlorocyclohexane isomers (other than $\gamma$ ), mg/kg, <i>Max</i>	0.3	
bb)	Lindane ( $\gamma$ -hexachlorocyclohexane), mg/kg, <i>Max</i>	0.6	
cc)	Malathion, mg/kg, <i>Max</i>	1.0	
dd)	Methidathion, mg/kg, <i>Max</i>	0.2	
ee)	Parathion, mg/kg, <i>Max</i>	0.5	
ff)	Parathion-methyl, mg/kg, <i>Max</i>	0.2	
gg)	Permethrin, mg/kg, <i>Max</i>	1.0	
hh)	Phosalone, mg/kg, <i>Max</i>	0.1	
jj)	Piperonyl butoxide, mg/kg, <i>Max</i>	3.0	
kk)	Pirimiphos-methyl, mg/kg, <i>Max</i>	4.0	
mm)	Pyrethrins (sum of), mg/kg, <i>Max</i>	3.0	
	Quintozene (sum of quintozene, pentachloroaniline and methyl pentachlorophenyl sulphide), mg/kg, <i>Max</i>	1.0	

## ANNEX A

(Clause 2)

## LIST OF REFERRED STANDARDS

<i>IS No.</i>	<i>Title</i>	<i>IS No.</i>	<i>Title</i>
IS 1070 : 2023	Reagent grade water — Specification ( <i>fourth revision</i> )		G <sub>1</sub> and G <sub>2</sub> in cereals, nuts and derived products — High performance liquid chromatographic method
IS 1797 : 2017	Spices and condiments — Methods of test ( <i>third revision</i> )	IS 16913 : 2018	Methods of test for cosmetics — Determination of heavy metals (arsenic, cadmium, lead and mercury) by atomic absorption spectrometry (AAS)
IS 4333 (Part 1) : 2018	Methods of analysis for foodgrains: Part 1 Refractions ( <i>third revision</i> )		
IS 13145 : 2014	Spices and condiments — Methods of sampling ( <i>second revision</i> )	IS 17924 : 2022	Determination of pesticide residue in herbal materials
IS 16287 : 2015/ ISO 16050 : 2003	Foodstuffs — Determination of aflatoxin B <sub>1</sub> , and the total content of aflatoxins B <sub>1</sub> , B <sub>2</sub> ,		

To access Indian Standards click on the link below:

[https://www.services.bis.gov.in/php/BIS\\_2.0/bisconnect/knownyourstandards/Indian\\_standards/isdetails/](https://www.services.bis.gov.in/php/BIS_2.0/bisconnect/knownyourstandards/Indian_standards/isdetails/)



## ANNEX B

[Table 1, SI No. vii), Col (4)]

## THIN-LAYER CHROMATOGRAM IDENTIFICATION

**B-1 GENERAL**

The chromatographic profile of the test solution is similar to that of the reference solution, prepared with the botanical reference standard BRS *Anjūr dashṭī* (fruit) under the same conditions, with respect to the position and fluorescence/colour of the bands.

**B-2 THIN LAYER CHROMATOGRAPHY (TLC)****B-2.1 Apparatus**

**B-2.1.1** *Thin Layer Chromatograph with a Pre Coated Silica Gel 60 F<sub>254</sub> Plate*

**B-2.1.2** *Analytical Balance*

**B-2.1.3** *Standard Glassware*

**B-2.1.4** *Water Bath*

**B-2.2 Reagents**

**B-2.2.1** *Preparation of Test Solution*

**B-2.2.2** *Reference Standard* — botanical reference standard

**B-2.2.3** *Chloroform* — AR or equivalent grade

**B-2.2.4** *Ethanol* — AR or equivalent grade

**B-2.2.5** *Toluene* — AR or equivalent grade

**B-2.2.6** *Ethyl Acetate* — AR or equivalent grade

**B-2.2.7** *Formic Acid* — AR or equivalent grade

**B-2.2.8** *Vanillin-Sulphuric Acid Reagent* — AR or equivalent grade

**B-2.3 Procedure****B-2.3.1 Solvent System**

Dissolve chloroform and toluene, ethyl acetate and formic acid in the ratio (8 : 2 : 0.1).

**B-2.3.2** Preparation of test solution to 3 g of the substance being examined, add 25 ml of methanol, heat on a water bath for 10 min to 15 min, cool and filter.

**B-2.2.3** Preparation of botanical reference standard (BRS) solution.

**B-3 IDENTIFICATION****B-3.1 Thin Layer Chromatography of Chloroform Extract**

Carry out thin layer chromatography on a pre-coated silica gel 60 F<sub>254</sub> plate. Test solution prepare the extract from 2 g of coarsely powdered plant material by refluxing with 20 ml of chloroform for 30 min. Filter, concentrate the extract and use as test solution.

**B-3.2 Procedure**

Apply 8 µl of the test solution as 8 mm band at a height of 10 mm from the base of a 10 cm × 10 cm TLC plate. Develop the plate to a distance of 8 cm from the band using solvent system – toluene : ethyl acetate : formic acid (8 : 2 : 0.1). Air dry the plate and examine under UV 254 nm and UV 366 nm. Spray the plate with vanillin-sulphuric acid reagent and heat at temperature of 105 °C to 110 °C till the bands appear. The chromatographic profile of the test solution is similar to that of the reference solution, prepared with the BRS *Anjūr dashṭī* (fruit) under the same conditions, with respect to the position and fluorescence/colour of the bands (see Fig. 2).

**B-3.3 Thin Layer Chromatography of Alcohol Extract**

Carry out thin layer chromatography on a pre-coated silica gel 60 F<sub>254</sub> plate.

**B-3.4 Test Solution**

Prepare the extract from 2 g of coarsely powdered plant material by refluxing with 20 ml of ethanol for 30 min. Filter, concentrate the extract and use as test solution.

**B-3.5 Procedure**

Apply 1 µl of the test solution as 8 mm band at a height of 10 mm from the base of a 10 cm × 10 cm TLC plate. Develop the plate to a distance of 8 cm from the band using solvent system – toluene : ethyl acetate : formic acid (7.8 : 2.2 : 0.1). Air dry the plate and examine under UV 366 nm. Spray the

plate with vanillin-sulphuric acid reagent and heat at temperature of 105 °C to 110 °C till the bands appear. The chromatographic profile of the test solution is similar to that of the reference solution, prepared with the BRS *Anjir dashi* (fruit) under

the same conditions, with respect to the position and fluorescence/colour of the bands (see Fig. 3).

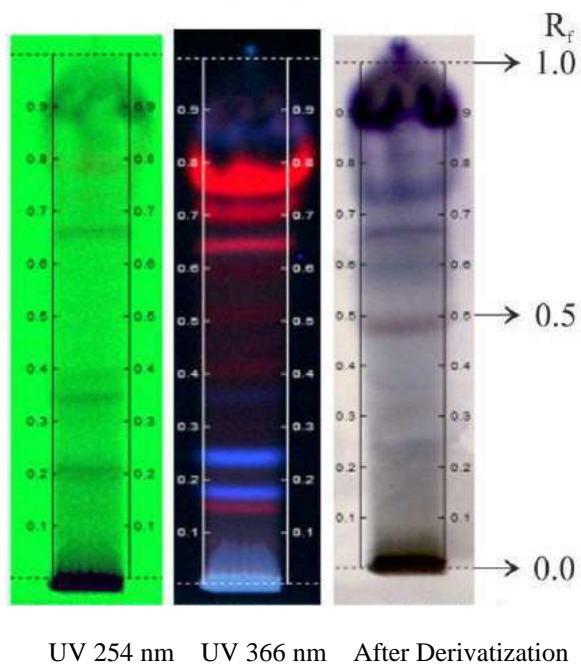


FIG. 2 TLC FINGERPRINTS OF ANJIR DASHTI (*Ficus hispida* L. f.) FRUIT

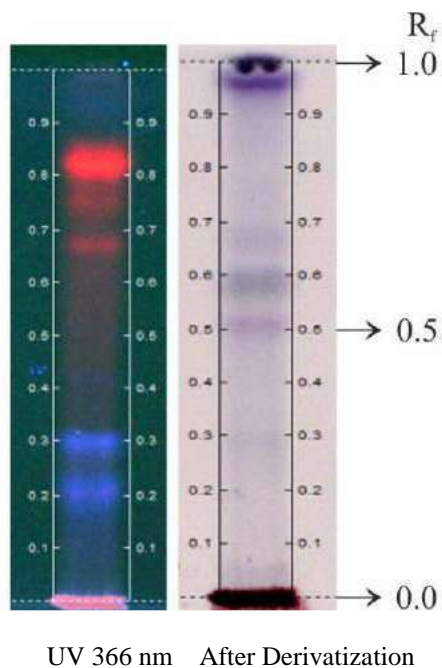


FIG. 3 TLC FINGERPRINTS OF ANJIR DASHTI (*Ficus hispida* L. f.) FRUIT

## ANNEX C

[Table 1, Sl No.(xiii), (xiv), (xv), (xvi), (xvii) and (xviii)]

## TEST FOR MICROBIAL LIMITS

**C-1 GENERAL**

The tests are designed for the estimation of the number of viable aerobic microorganisms present and for detecting the presence of designated microbial species in the extract. The term 'growth' is used to designate the presence and presumed proliferation of viable microorganisms.

**C-2 APPARATUS****C-2.1 Oven for Dry Sterilization****C-2.2 Autoclave for Wet Sterilization****C-2.3 Incubator****C-2.4 Water Bath****C-2.5 pH-Meter**

**C-2.6 Sterile Membrane Filters** — 50 mm in diameter

**C-2.7 Colony Counting Equipment****C-2.8 Analytical Balance****C-2.9 Standard Glassware****C-3 MEDIA**

Culture media may be prepared as given below or dehydrated culture media may be used provided that, when reconstituted as directed by the manufacturer, they have similar ingredients and/or yield media comparable to those obtained from the formulae given below.

Where agar is specified in a formula, use agar that has a moisture content of not more than 15 percent. Where water is called for in a formula, use purified water. Unless otherwise indicated, the media should be sterilized by heating in an autoclave (15 psi) at 121 °C for 15 min. In preparing media by the formulas given below, dissolve the soluble solids in the water, using heat if necessary, to effect complete solution, add solutions of 0.1 N hydrochloric acid or 0.1N sodium hydroxide in quantities sufficient to yield the required pH in the medium when it is ready for use. Determine the pH at 25 °C ± 2 °C.

**C-3.1 Baird Parker Agar Medium****C-3.1.1 Composition**

Pancreatic digest of casein	10.0 g
Beef extract	5.0 g
Yeast extract	1.0 g
Lithium chloride	5.0 g
Agar	20.0 g
Glycine	12.0 g
Sodium pyruvate	10.0 g
Water	1 000 ml

**C-3.1.2 Preparation**

Suspend the components in 1 000 ml of water, heat with frequent agitation and boil for 1 min. Sterilize, cool in between 45 °C to 50 °C, add 10 ml of a one percent w/v solution of sterile potassium tellurite and 50 ml of egg yolk emulsion. Mix thoroughly, but gently and pour into plates. (Prepare the egg yolk emulsion by disinfecting the surface of whole shell eggs, aseptically cracking the eggs, and separating out intact yolks into a sterile graduated cylinder. Add sterile saline solution, get a 3 to 7 ratio of egg yolk to saline. Add to a sterile blender cup and mix at high speed for 5 seconds). Adjust the pH after sterilization to 6.8 ± 0.2.

**C-3.2 Bismuth Sulphite Agar Medium****C-3.2.1 Composition**

Solution (1)	
Beef extract	6 g
Peptone	10 g
Agar	24 g
Ferric citrate	0.4 g
Brilliant green	10 mg
Water	1 000 ml
Solution (2)	
Ammonium bismuth citrate	3 g
Sodium sulphite	10 g
Anhydrous disodium hydrogen phosphate	5 g
Dextrose monohydrate	5 g
Water	100 ml

**C-3.2.2 Preparation**

Suspend the components of Solution (1) in 1 000 ml of water. Heat to boiling to dissolve the medium completely. Sterilize by maintaining at 115 °C for 30 min.

Suspend the components of Solution (2) in 100 ml of water. Heat to boiling to dissolve the medium completely. Do not autoclave.

Add 1 volume of solution (2) to 10 volumes of solution (1) previously melted and cooled to a temperature of 55 °C. Bismuth sulphite agar medium should be stored at 2 °C to 8 °C for 5 days before use.

**C-3.3 Brilliant Green Agar Medium****C-3.3.1 Composition**

Peptone	10.0 g
Yeast extract	3.0 g
Lactose	10.0 g
Sucrose	10.0 g
Sodium chloride	5.0 g
Phenol red	80.0 g
Brilliant green	12.5 mg
Agar	12.0 g
Water	1 000 ml

**C-3.3.2 Preparation**

Mix the components, allow to stand for 15 min, sterilize by maintaining at 115 °C for 30 min and mix before pouring.

**C-3.4 Buffered Sodium Chloride Peptone Solution pH 7.0****C-3.4.1 Composition**

Potassium dihydrogen phosphate	3.56 g
Disodium hydrogen phosphate	7.23 g
Sodium chloride	4.30 g
Peptone (meat or casein)	1.0 g
Water	1 000 ml

**C-3.4.2 Preparation**

Mix the components and heat if necessary to dissolve the medium completely. 0.1 percent to 1.0 percent *w/v* polysorbate 20 or polysorbate 80 may be added. Sterilize by heating in an autoclave at 121 °C for 15 min.

**C-3.5 Casein Soyabean Digest Agar Medium****C-3.5.1 Composition**

Pancreatic digest of casein	15.0 g
Papaic digest of soyabean meal	5.0 g
Sodium chloride	5.0 g
Agar	15.0 g
Water	1 000 ml

**C-3.5.2 Preparation**

Suspend the components in water. Heat to boiling to dissolve the medium completely. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to 7.3 ± 0.2.

**C-3.6 Cetrinide Agar Medium****C-3.6.1 Composition**

Pancreatic digest of gelatin	20.0 g
Magnesium chloride	1.4 g
Potassium sulphate	10.0 g
Cetrinide	0.3 g
Agar	13.6 g
Glycerin	10.0 g
Water	1 000 ml

**C-3.6.2 Preparation**

Dissolve the components in 1 000 ml of water. Heat to boiling for 1 min with shaking. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to 7.0 to 7.4.

**C-3.7 Deoxycholate Citrate Agar Medium****C-3.7.1 Composition**

Beef extract	5.0 g
Peptone	5.0 g
Lactose	10.0 g
Trisodium citrate	8.5 g
Sodium thiosulphate	5.4 g
Ferric citrate	1.0 g
Sodium deoxycholate	5.0 g
Neutral red	0.02 g
Agar	12.0 g
Water	1 000 ml

**C-3.7.2 Preparation**

Mix the components and allow to stand for 15 min. Gently boil with continuous stirring and continue boiling until solution is complete. Cool to 80 °C, mix, pour and cool rapidly.

Care should be taken not to overheat deoxycholate citrate agar during preparation. It should not be re-melted and the surface of the plates should be dried before use.

**C-3.8 Fluid Casein Digest Soya Lecithin Polysorbate 20 Medium****C-3.8.1 Composition**

Pancreatic digest of casein	20 g
Soya lecithin	5 g
Polysorbate 20	40 ml
Water	1 000 ml

**C-3.8.2 Preparation**

Dissolve the pancreatic digest of casein and soya lecithin in water, heating in a water bath at 48 °C to 50 °C for about 30 min to effect solution. Add polysorbate 20, mix and dispense as desired. Sterilize at 121 °C for 15 min in an autoclave.

**C-3.9 Fluid Lactose Medium****C-3.9.1 Composition**

Beef extract	3.0 g
Pancreatic digest of gelatin	5.0 g
Lactose	5.0 g
Water	1 000 ml

**C-3.9.2 Preparation**

Suspend the components in 1 000 ml water. Heat if necessary to dissolve the medium completely. Sterilize at 121 °C for 15 min in an autoclave. Cool as quickly as possible after sterilization. Adjust the pH after sterilization to  $6.9 \pm 0.2$ .

**C-3.10 Lactose Broth Medium****C-3.10.1 Composition**

Beef extract	3.0 g
Pancreatic digest of gelatin	5.0 g
Lactose	5.0 g
Water	1 000 ml

**C-3.10.2 Preparation**

Suspend the components in water and heat if necessary to dissolve the medium completely. Sterilize at 121 °C for 15 min in an autoclave. Adjust the pH after sterilization to  $6.9 \pm 0.2$ .

**C-3.11 Levine Eosin Methylene Blue Agar Medium****C-3.11.1 Composition**

Pancreatic digest of gelatin	10.0 g
Dibasic potassium phosphate	2.0 g
Agar	15.0 g
Lactose	10.0 g
Eosin Y	400 mg
Methylene blue	65 mg
Water	1 000 ml

**C-3.11.2 Preparation**

Dissolve the pancreatic digest of gelatin, dibasic potassium phosphate and agar in water with warming and allow to cool. Just prior to use, liquify the gelled agar solution and the remaining ingredients, as solutions, in the following amounts and mix. For each 100 ml of the liquified agar solution use 5 ml of a 20 percent w/v solution of lactose, 2 ml of a 2 percent w/v solution of eosin Y and 2 ml of a 0.33 percent w/v solution of methylene blue. The finished medium may not be clear. Adjust the pH after sterilization to  $7.1 \pm 0.2$ .

**C-3.12 MacConkey Agar Medium****C-3.12.1 Composition**

Pancreatic digest of gelatin	17.0 g
Peptone (meat and casein, equal parts)	3.0 g
Lactose	10.0 g
Sodium chloride	5.0 g
Bile salts	1.5 g
Agar	13.5 g
Neutral red	30 mg
Crystal violet	1 mg
Water	1 000 ml

**C-3.12.2 Preparation**

Boil the mixture of solids and water for 1 min to effect solution. Sterilize at 121 °C for 15 min in an autoclave. Adjust the pH after sterilization to  $7.1 \pm 0.2$ .

**C-3.13 MacConkey Broth Medium****C-3.13.1 Composition**

Pancreatic digest of gelatin	20.0 g
Lactose	10.0 g
Dehydrated ox bile	5.0 g
Bromocresol purple	10 mg
Water	1 000 ml

**C-3.13.2 Preparation**

Suspend the components in 1 000 ml of water and heat if necessary to dissolve the medium completely. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to  $7.3 \pm 0.2$ .

**C-3.14 Mannitol Salt Agar Medium****C-3.14.1 Composition**

Pancreatic digest of gelatin	5.0 g
Peptic digest of animal tissue	5.0 g
Beef extract	1.0 g
D-Mannitol	10.0 g
Sodium chloride	75.0 g
Agar	15.0 g
Phenol red	25 mg
Water	1 000 ml

**C-3.14.2 Preparation**

Mix the components, heat with frequent agitation and boil for 1 min to effect solution. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to  $7.4 \pm 0.2$ .

**C-3.15 Nutrient Broth Medium****C-3.15.1 Composition**

Beef extract	10.0 g
Peptone	10.0 g
Sodium chloride	5 mg
Water	1 000 ml

**C-3.15.2 Preparation**

Dissolve the components with the aid of heat. Adjust the pH to 8.0 to 8.4 with 5 M sodium hydroxide and boil for 10 min. Filter and sterilize by maintaining at 115 °C for 30 min and adjust pH to  $7.3 \pm 0.1$ .

**C-3.16 Nutrient Agar Medium**

Nutrient broth gelled by the addition of 1 percent to 2 percent w/v of agar.

**C-3.17 Pseudomonas Agar Medium for Detection of Fluorescein****C-3.17.1 Composition**

Pancreatic digest of casein	10.0 g
Peptic digest of animal tissue	10.0 g
Anhydrous dibasic potassium phosphate	1.5 g
Magnesium sulphate hepta hydrate	1.5 g
Glycerin	10.0 ml
Agar	15.0 g
Water	1 000 ml

**C-3.17.2 Preparation**

Dissolve the solid components in water before adding glycerin. Heat with frequent agitation and boil for 1 min to effect solution. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to  $7.2 \pm 0.2$ .

**C-3.18 Pseudomonas Agar Medium for Detection of Pyocyanin****C-3.18.1 Composition**

Pancreatic digest of gelatin	20.0 g
Anhydrous magnesium chloride	1.4 g
Anhydrous potassium sulphate	10.0 g
Agar	15.0 g
Glycerin	10.0 ml
Water	1 000 ml

**C-3.18.2 Preparation**

Dissolve the solid components in water before adding glycerin. Heat with frequent agitation and boil for 1 min to effect solution. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to  $7.2 \pm 0.2$ .

**C-3.19 Sabouraud Dextrose Agar Medium****C-3.19.1 Composition**

Dextrose	40 g
Peptic digest of animal tissue and pancreatic digest of casein (1 : 1)	10 g
Agar	15 g
Water	1 000 ml

**C-3.19.2 Preparation**

Mix the components and heat to boiling to dissolve completely. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to  $5.6 \pm 0.2$ .

**C-3.20 Sabouraud Dextrose Agar Medium with Antibiotics**

To 1 000 ml of Sabouraud dextrose agar medium, add 0.1 g of benzylpenicillin sodium and 0.1 g of tetracycline HCl or alternatively add 50 mg of chloramphenicol immediately before use.

**C-3.21 Selenite F Broth****C-3.21.1 Composition**

Peptone	5 g
Lactose	4 g
Disodium hydrogen phosphate	10 g
Sodium hydrogen selenite	4 g
Water	1 000 ml

**C-3.21.2 Preparation**

Suspend the components in water and mix well. Warm to dissolve the medium completely. Distribute in sterile containers and sterilize by maintaining at 100 °C for 30 min.

**C-3.22 Fluid Selenite Cystine Medium****C-3.22.1 Composition**

Pancreatic digest of casein	5.0 g
Lactose	4.0 g
Sodium phosphate	10.0 g
Sodium hydrogen selenite	4.0 g
l-Cystine	10.0 mg
Water	1 000 ml

**C-3.22.2 Preparation**

Suspend the components in water and heat in flowing steam for 15 min. Adjust the final pH to 7.0 ± 0.2. Do not sterilize.

**C-3.23 Tetrathionate Broth Medium****C-3.23.1 Composition**

Beef extract	0.9 g
Peptone	4.5 g
Yeast extract	1.8 g
Sodium chloride	4.5 g
Calcium carbonate	25.0 g
Sodium thiosulphate	40.7 g
Water	1 000 ml

**C-3.23.2 Preparation**

Dissolve the solids in 1 000 ml of water and heat the solution to boil. On the day of use, add a

solution prepared by dissolving 5 g of potassium iodide and 6 g of iodine in 20 ml of water.

**C-3.24 Tetrathionate Bile Brilliant Green Broth Medium****C-3.24.1 Composition**

Peptone	8.6 g
Dehydrated ox bile	8.0 g
Sodium chloride	6.4 g
Calcium carbonate	20.0 g
Potassium tetrathionate	20.0 g
Brilliant green	70 mg
Water	1 000 ml

**C-3.24.2 Preparation**

Suspend the components in 1 000 ml of water. Heat just to boiling. Do not autoclave or reheat. Adjust the pH so that after heating it is 7.0 ± 0.2.

**C-3.25 Triple Sugar Iron Agar Medium****C-3.25.1 Composition**

Beef extract	3.0 g
Yeast extract	3.0 g
Peptone	20.0 g
Lactose	10.0 g
Sucrose	10.0 g
Dextrose monohydrate	1.0 g
Ferrous sulphate	0.2 g
Sodium chloride	5.0 g
Sodium thiosulphate	0.3 g
Phenol red	24 mg
Water	1 000 ml

**C-3.25.2 Preparation**

Mix the components in 1 000 ml of water, allow standing for 15 min. Heat to boiling to dissolve the medium completely. Mix well and distribute into test tubes and Sterilize by maintaining at 121 °C for 15 min. Allow the medium to stand in a sloped form with a butt about 2.5 cm long.

**C-3.26 Urea Broth Medium****C-3.26.1 Composition**

Potassium dihydrogen orthophosphate	9.1 g
Anhydrous disodium hydrogen phosphate	9.5 g
Urea	20.0 g
Yeast extract	0.1 g

## IS 18782 : 2024

Phenol red	10 mg
Water	1 000 ml

### C-3.26.2 Preparation

Mix the components, sterilize by filtration and distribute aseptically in sterile containers.

### C-3.27 Vogel Johnson Agar Medium

#### C-3.27.1 Composition

Pancreatic digest of casein	10.0 g
Yeast extract	5.0 g
Mannitol	10.0 g
Dibasic potassium phosphate	5.0 g
Lithium chloride	5.0 g
Glycerin	10.0 g
Agar	16.0 g
Phenol red	25.0 mg
Water	1 000 ml

#### C-3.27.2 Preparation

Suspend the components in 1 000 ml of water. Boil the solution of solids for 1 min. Sterilize at 121 °C for 15 min in an autoclave. Cool to 45 °C to 50 °C and add 20 ml of 1 percent w/v sterile solution of potassium tellurite. Adjust the pH after sterilization to 7.0 ± 0.2.

### C-3.28 Xylose Lysine Deoxycholate Agar Medium

#### C-3.28.1 Composition

Xylose	3.5 g
l-Lysine	5.0 g
Lactose	7.5 g
Sucrose	7.5 g
Sodium chloride	5.0 g
Yeast extract	3.0 g
Phenol red	80 mg
Agar	13.5 g
Sodium deoxycholate	2.5 g
Sodium thiosulphate	6.8 g
Ferric ammonium citrate	800 mg
Water	1 000 ml

#### C-3.28.2 Preparation

Suspend the components in 1 000 ml of water. Heat with frequent agitation until the medium boils. Do not overheat or sterilize. Transfer at once to a water bath maintained at about 50 °C and pour into plates

as soon as the medium has cooled. Adjust the final pH to 7.4 ± 0.2.

### C-4 SAMPLING

Use 10 ml or 10 g specimens for each of the tests specified in this Indian Standard.

### C-5 PRECAUTIONS

The microbial limit tests should be carried out under conditions designed to avoid accidental contamination during the test. The precautions taken to avoid contamination must be such that they do not adversely affect any microorganisms that should be revealed in the test.

### C-6 PRELIMINARY TESTING

Prior to doing the tests, inoculate diluted extracts being examined with separate viable cultures of *Escherichia coli*, *Salmonella* species, *Pseudomonas aeruginosa* and *Staphylococcus aureus*. This is done by adding 1 ml of 24 h broth culture containing not less than 1 000 microorganisms to the first dilution (in buffer solution pH 7.2, fluid soya bean casein digest medium or fluid lactose medium) of the test material and following the test procedure. If the organisms fail to grow in the relevant medium the procedure should be modified by (a) increasing the volume of diluent with the quantity of test material remaining the same, or (b) incorporating a sufficient quantity of a suitable inactivating agent in the diluents, or (c) combining the afore mentioned modifications so as to permit growth of the organisms in the media. If inhibitory substances are present in the extracts, 0.5 percent of soya lecithin and 4 percent of polysorbate 20 may be added to the culture medium. Alternatively, repeat the test as described in the previous paragraph, using fluid casein digest soya lecithin polysorbate 20 medium, to demonstrate neutralization of preservatives or other antimicrobial agents in the test material. Where inhibitory substances are contained in the extracts and the latter is soluble, the membrane filtration method described under total aerobic microbial count may be used. If in spite of incorporation of suitable inactivating agents and a substantial increase in the volume of diluent it is still not possible to recover the viable cultures described above and where the article is not suitable for applying the membrane filtration method, it can be assumed that the failure to isolate the inoculated organism may be due to the bactericidal activity of the product. This may indicate that the article is not likely to be contaminated with the given species of microorganisms. However, monitoring should be continued to establish the spectrum of inhibition and bactericidal activity of the article.



## C-7 TOTAL AEROBIC MICROBIAL COUNT

Pre-treat the extracts and raw materials being examined as described below.

NOTE — The raw materials need to be ground as a coarse powder before analysis.

### C-7.1 Water Soluble Products

Dissolve 10 g or dilute 10 ml of the extract preparation being examined, unless otherwise specified, in buffered sodium chloride peptone solution pH 7.0 or any other suitable medium shown to have no antimicrobial activity under the conditions of test and adjust the volume to 100 ml with the same medium. If necessary, adjust the pH to about 7.

### C-7.2 Products Insoluble in Water (Non Fatty)

Suspend 10 g or 10 ml of the extract preparation being examined, unless otherwise specified, in buffered sodium chloride peptone solution pH 7.0 or any other suitable medium shown not to have antimicrobial activity under the conditions of the test and dilute to 100 ml with the same medium. If necessary, divide the preparation being examined and homogenize the suspension mechanically. A suitable surface active agent such as 0.1 percent w/v of polysorbate 80 may be added to assist the suspension of poorly wettable substances. If necessary, adjust pH of the suspension to about 7.

### C-7.3 Fatty Products

Homogenise 10 g or 10 ml of the extract preparation being examined, unless otherwise specified, with 5 g of polysorbate 20 or polysorbate 80. If necessary, heat to not more than 40 °C. Mix carefully while maintaining the temperature in the water bath or in an oven. Add 85 ml of buffered sodium chloride peptone solution pH 7.0 or any other suitable medium shown to have no antimicrobial activity under the conditions of the test, heated to not more than 40 °C if necessary. Maintain this temperature for the shortest time necessary for formation of an emulsion and in any case for not more than 30 min. If necessary, adjust the pH to about 7.

### C-7.4 Examination of the Sample

Determine the total aerobic microbial count in the extract being examined by any of the following methods.

### C-7.5 Membrane Filtration

Use membrane filters 50 mm in diameter and having a nominal pore size not greater than

0.45 µm the effectiveness of which in retaining bacteria has been established for the type of preparation being examined. Transfer 10 ml or a quantity of each dilution containing 1g of the preparation being examined to each of two membrane filters and filter immediately. If necessary, dilute the pretreated extract preparation so that a colony count of 10 to 100 may be expected. Wash each membrane by filtering through it three or more successive quantities, each of about 100 ml, of a suitable liquid such as buffered sodium chloride peptone solution pH 7.0. For fatty substances add to the liquid polysorbate 20 or polysorbate 80. Transfer one of the membrane filters, intended for the enumeration of bacteria, to the surface of a plate of casein soyabean digest agar and the other, intended for the enumeration of fungi, to the surface of a plate of Sabouraud dextrose agar with antibiotics.

Incubate the plates for 5 days, unless a more reliable count is obtained in shorter time, at 30 °C to 35 °C in the test for bacteria and 20 °C to 25 °C in the test for fungi. Count the number of colonies that are formed. Calculate the number of microorganisms per g or per ml of the extract preparation being examined, if necessary count bacteria and fungi separately.

### C-7.6 Plate Count for Bacteria

Using Petri dishes 9 to 10 cm in diameter, add to each dish a mixture of 1 ml of the pretreated extract preparation and about 15 ml of liquified casein soyabean digest agar at not more than 45 °C.

Alternatively, spread the pretreated extract preparation on the surface of the solidified medium in a petri dish of the same diameter. If necessary, dilute the pretreated extract preparation as described above so that a colony count of not more than 300 may be expected. Prepare at least two such petri dishes using the same dilution and incubate at 30 °C to 35 °C for 5 days, unless a more reliable count is obtained in a shorter time. Count the number of colonies that are formed. Calculate the results using plates with the greatest number of colonies but taking 300 colonies per plate as the maximum consistent with good evaluation.

### C-7.7 Plate Count for Fungi

Proceed as described in the test for bacteria but use Sabouraud dextrose agar with antibiotics in place of casein soyabean digest agar and incubate the plates at 20 °C to 25 °C for 5 days, unless a more reliable count is obtained in a shorter time. Calculate the results using plates with not more than 100 colonies.

**C-7.8 Multiple Tube or Serial Dilution Method**

In each of fourteen test tubes of similar size place 9.0 ml of sterile fluid soyabean casein digest medium. Arrange twelve of the tubes in four sets of three tubes each. Put aside one set of three tubes to serve as controls. Into each of three tubes of one set ("100") and into fourth tube (A) pipette 1 ml of the solution of suspension of the test specimen (extract) and mix. From tube A pipette 1 ml of its contents into the one remaining tube (B) not included in the set and mix. These two tubes contain 100 mg (or 100 µl) and 10 mg (or 10 µl) of the specimen respectively. Into each of the second set ("10") of three tubes pipette 1 ml from tube A, and into each tube of the third set ("1") pipette 1 ml from tube B. Discard the unused contents of tube A and tube B. Close well and incubate all of the tubes. Following the incubation period, examine the tubes for growth. The three control tubes remain clear. Observations in the tubes containing the test specimen, when interpreted by reference to [Table 2](#), indicate the most probable number of

microorganisms per g or per ml of the test specimen.

**C-8 TESTS FOR SPECIFIED MICROORGANISMS****C-8.1 Pre-Treatment of the Extract Sample Being Examined**

Proceed as described under the test for total aerobic microbial count but using lactose broth or any other suitable medium shown to have no antimicrobial activity under the conditions of test in place of buffered sodium chloride peptone solution pH 7.0.

**C-8.2 *Escherichia coli***

Place the prescribed quantity in a sterile screw capped container, add 50 ml of nutrient broth, shake, allow to stand for 1 h (4 h for gelatin) and shake again. Loosen the cap and incubate at 37 °C for 18 h to 24 h.

**Table 2 Most probable Total Count by Multiple Tube or Serial Dilution Method**(Clause [C-7.8](#))

Sl No.	Observed Combination of Numbers of Tubes Showing Growth in Each Set			Most Probable Number of Microorganisms per g or per ml
	Number of mg (or ml) of Specimen per Tube			
	100 (100 µl)	10 (10 µl)	1 (1 µl)	
(1)	(2)	(3)	(4)	(5)
i)	3	3	3	> 1 100
ii)	3	3	2	1 100
iii)	3	3	1	500
iv)	3	3	0	200
v)	3	2	3	290
vi)	3	2	2	210
vii)	3	2	1	150
viii)	3	2	0	90
ix)	3	1	3	160
x)	3	1	2	120
xi)	3	1	1	70
xii)	3	1	0	40
xiii)	3	0	3	95
xiv)	3	0	2	60
xv)	3	0	1	40
xvi)	3	0	0	23

**C-8.2.1 Primary Test**

Add 1.0 ml of the enrichment culture to a tube containing 5 ml of MacConkey broth. Incubate in a water bath at 36 °C to 38 °C for 48 h. If the contents of the tube show acid and gas, carry out the secondary test.

**C-8.2.2 Secondary Test**

Add 0.1 ml of the contents of the tubes containing (a) 5 ml of MacConkey broth and (b) 5 ml of peptone water. Incubate in a water bath at 43.5 °C to 44.5 °C for 24 h and examine tube (a) for acid and gas and tube (b) for indole. To test for indole, add 0.5 ml of Kovac's reagent, shake well and allow to stand for 1 min. If a red colour is produced in the reagent layer indole is present. The presence of acid and gas and of indole in the secondary test indicates the presence of *Escherichia coli*.

Carry out a control test by repeating the primary and secondary tests, adding 1.0 ml of the enrichment culture and a volume of broth containing 10 to 50 *Escherichia coli* (NCTC 9002) organisms, prepared from a 24 h culture in nutrient broth, to 5 ml of MacConkey broth. The test is not valid unless the results indicate that the control contains *Escherichia coli*.

**C-8.2.3 Alternative Test**

By means of an inoculating loop, streak a portion from the enrichment culture (obtained in the previous test) on the surface of MacConkey agar medium. Cover and invert the dishes and incubate.

Upon examination, if none of the colonies are brick red in colour and have a surrounding zone of precipitated bile the sample meets the requirements

of the test for the absence of *Escherichia coli*. If the colonies described above are found, transfer the suspect colonies individually to the surface of Levine eosin methylene blue agar medium, plated on petri dishes. Cover and invert the plates and incubate. Upon examination, if none of the colonies exhibits both a characteristic metallic sheen under reflected light and a blue black appearance under transmitted light, the sample meets the requirements of the test for the absence of *Escherichia coli*. The presence of *Escherichia coli* may be confirmed by further suitable cultural and biochemical tests.

**C-8.3 Salmonella**

Transfer a quantity of the pretreated preparation being examined containing 1 g or 1 ml of the product to 100 ml of nutrient broth in a sterile screw capped jar, shake, allow to stand for 4 h and shake again. Loosen the cap and incubate at 35 °C to 37 °C for 24 h.

**C-8.3.1 Primary Test**

Add 1.0 ml of the enrichment culture to each of the two tubes containing (a) 10 ml of selenite F broth and (b) tetrathionate bile brilliant green broth and incubate at 36 °C to 38 °C for 48 h. From each of these two cultures subculture on at least two of the following four agar media: bismuth sulphate agar, brilliant green agar, deoxycholate citrate agar and xylose lysine deoxycholate agar. Incubate the plates at 36 °C to 38 °C for 18 to 24 h. Upon examination, if none of the colonies conforms to the description given in [Table 3](#), the sample meets the requirements of the test for the absence of the genus *Salmonella*. If any colonies conforming to the description in [Table 3](#) are produced, carry out the secondary test.

**Table 3 Interpretation of the Confirmatory Tests for *Salmonella***(Clauses [C-8.3.1](#) and [C-8.3.2](#))

Sl No.	Medium	Description of Colony
(1)	(2)	(3)
i)	Bismuth sulphite agar	Black or green
ii)	Brilliant green agar	Small, transparent and colorless, or opaque, pinkish or white (frequently surrounded by a pink or red zone)
iii)	Deoxycholate citrate agar	Colorless and opaque, with or without black centers
iv)	Xylose lysine deoxycholate agar	Red with or without black centres

**C-8.3.2 Secondary Test**

Subculture any colonies showing the characteristics given in [Table 3](#) in triple sugar iron agar by first inoculating the surface of the slope and then making a stab culture with the same inoculating needle, and at the same time inoculate a tube of urea broth. Incubate at 36 °C to 38 °C for 18 h to 24 h. The formation of acid and gas in the stab culture (with or without concomitant blackening) and the absence of acidity from the surface growth in the triple sugar iron agar, together with the absence of a red color in the urea broth, indicate the presence of *Salmonella*. If acid but no gas is produced in the stab culture, the identity of the organisms should be confirmed by agglutination tests.

Carry out the control test by repeating the primary and secondary tests using 1.0 ml of the enrichment culture and a volume of broth containing 10 to 50 *Salmonella abony* (NCTC 6017) organisms, prepared from a 24 h culture in nutrient broth, for the inoculation of the tube (a) and tube (b). The test is not valid unless the results indicate that the control contains *Salmonella*.

**A-8.4 *Pseudomonas aeruginosa***

Pretreat the preparation being examined as described above and inoculate 100 ml of fluid soyabean casein digest medium with a quantity of the solution, suspension or emulsion thus obtained containing 1 g or 1 ml of the preparation being examined. Mix and incubate at 35 °C to 37 °C for 24 to 48 h. Examine the medium for growth and if growth is present, streak a portion of the medium on the surface of cetrimide agar medium, each plated on petri dishes. Cover and incubate at 35 °C to 37 °C for 18 h to 24 h.

If, upon examination, none of the plates contains colonies having the characteristics listed in [Table 4](#) for the media used, the sample meets the requirement for freedom from *Pseudomonas aeruginosa*. If any colonies conforming to the description in [Table 4](#) are produced, carry out the oxidase and pigment tests. Streak representative suspect colonies from the agar surface of cetrimide agar on the surfaces of *Pseudomonas* agar medium for detection of fluorescein and *Pseudomonas* agar medium for detection of pyocyanin contained in petri dishes. Cover and invert the inoculated media and incubate at 33 °C to 37 °C for not less than 3 days. Examine the streaked surfaces under ultra violet light. Examine the plates to determine whether colonies conforming to the description in [Table 4](#) are present. If growth of suspect colonies occurs, place 2 or 3 drops of a freshly prepared 1 percent *w/v* solution of N,N,N1,N1-tetramethyl-4-phenylenediamine dihydrochloride on filter paper and smear with the colony. If there is no development of a pink color, changing to purple, the sample meets the requirements of the test for the absence of *Pseudomonas aeruginosa*.

**C-8.5 *Staphylococcus aureus***

Proceed as described under *Pseudomonas aeruginosa*, if upon examination of the incubated plates, none of them contains colonies having the characteristics listed in for the media used, the sample meets the requirements for the absence of *Staphylococcus aureus*. If growth occurs, carry out the coagulase test. Transfer representative suspect colonies from the agar surface of any of the media listed in [Table 5](#) to individual tubes, each containing 0.5 ml of mammalian, preferably rabbit or horse, plasma with or without additives.

**Table 4 Interpretation of the Confirmatory Tests for *Pseudomonas aeruginosa***(Clause [C-8.4](#))

Sl No.	Selective Medium	Characteristic Colonial Morphology	Fluorescence in UV Light	Oxidase Test	Gram Stain
(1)	(2)	(3)	(4)	(5)	(6)
i)	Cetrimide agar	Generally greenish	Greenish	Positive	Negative rods
ii)	<i>Pseudomonas</i> agar medium for detection n of fluorescein	Generally colorless to yellowish	Yellowish	Positive	Negative rods
iii)	<i>Pseudomonas</i> agar for detection of pyocyanin	Generally greenish	Blue	Positive	Negative rods

**Table 5 Interpretation of the Confirmatory Tests for *Staphylococcus aureus***(Clause [C-8.5](#))

SI No.	Selective Medium	Characteristic Colonial Morphology	Gram Stain
(1)	(2)	(3)	(4)
i)	Vogel Johnson agar	Black surrounded by yellow zones	Positive cocci (in clusters)
ii)	Mannitol salt agar	Yellow colonies with yellow zones	Positive cocci (in clusters)
iii)	Baird parker agar	Black, shiny, surrounded by clear zones of 2 mm to 5 mm	Positive cocci (in clusters)

Incubate in water bath at 37 °C examining the tubes at 3 h and subsequently at suitable intervals up to 24 h. If no coagulation in any degree is observed, the sample meets the requirements of the test for the absence of *Staphylococcus aureus*.

#### **C-8.6 Validity of the Tests for Total Aerobic Microbial Count**

Grow the following test strains separately in tubes containing fluid soyabean casein digest medium at 30 °C to 35 °C for 18 to 24 h or, for *Candida albicans*, at 20 °C for 48 h.

*Staphylococcus aureus* (ATCC 6538; NCTC 10788)  
*Bacillus subtilis* (ATCC 6633; NCIB 8054)  
*Escherichia coli* (ATCC 8739; NCIB 8545)  
*Candida albicans* (ATCC 2091; ATCC 10231)

Dilute portions of each of the cultures using buffered sodium chloride peptone solution pH 7.0 to make test suspensions containing about 100 viable microorganisms per ml. Use the suspension of each of the microorganisms separately as a control of the counting methods, in the presence and absence of the preparation being examined, if necessary.

A count for any of the test organisms differing by not more than a factor of 10 from the calculated value for the inoculum should be obtained. To test the sterility of the medium and of the diluent and the aseptic performance of the test, carry out the total aerobic microbial count method using sterile buffered sodium chloride peptone solution pH 7.0 as the test preparation. There should be no growth of microorganisms.

#### **C-8.7 Validity of the Tests for Specified Microorganisms**

Grow separately the test strains of *Staphylococcus aureus* and *Pseudomonas aeruginosa* in fluid soyabean casein digest medium and *Escherichia coli* and *Salmonella typhimurium* at 30 °C to 35 °C for 18 h to 24 h. Dilute portions of each of the cultures using buffered sodium chloride peptone solution pH 7.0 to make test suspensions containing about 10<sup>3</sup> viable microorganisms per ml. Mix equal volume of each suspension and use 0.4 ml (approximately 10<sup>2</sup> microorganisms of each strain) as an inoculum in the test for *E. coli*, *S. typhimurium*, *P. aeruginosa* and *S. aureus*, in the presence and absence of the extract preparation being examined, if necessary. A positive result for the respective strain of microorganism should be obtained.

## ANNEX D

*(Foreword)*

## COMMITTEE COMPOSITION

Unani Sectional Committee, AYD 04

<i>Organization</i>	<i>Representative(s)</i>
School of Unani Medicine Education and Research, (SUMER), Jamia Hamdard, New Delhi	PROF ASIM ALI KHAN ( <i>Chairperson</i> )
Aligarh Muslim University, Aligarh	DR FAROOQ AHMAD DAR DR MOHAMMAD MOHSIN ( <i>Alternate I</i> ) DR SUMBUL REHMAN ( <i>Alternate II</i> )
Ayurvedic and Unani Tibbia College and Hospital, New Delhi	PROF MOHAMMED ZUBAIR DR MOHD FAROOQUE ( <i>Alternate I</i> ) DR MOHD NAUMAN SALEEM ( <i>Alternate II</i> )
Central Council for Research in Unani Medicine, New Delhi	DR RAM PRATAP MEENA DR FARAH AHMED ( <i>Alternate I</i> ) DR RITU KARWASRA ( <i>Alternate II</i> )
Central Government Health Scheme, Ministry of Health and Family Welfare, Govt of India, New Delhi	DR ABDUL QAYYUM DR MUZAMIL REHMAN ( <i>Alternate I</i> ) DR SUHAIL AKHTAR ( <i>Alternate II</i> )
Delhi Pharmaceutical Sciences and Research University, New Delhi	PROF AJAY SHARMA DR MUKESH NANDAVE ( <i>Alternate I</i> ) DR ARYA LAKSHMI MARISSETTI ( <i>Alternate II</i> )
Directorate of Ayush, Govt of NCT, New Delhi	DR SHAGUFA NASREEN DR PARAS WANI ( <i>Alternate I</i> ) DR FARAH NAAZ ( <i>Alternate II</i> )
Hamdard Laboratories, New Delhi	DR SANTOSH KUMAR JOSHI DR SAGHEER AHMAD KHAN ( <i>Alternate</i> )
Himalaya Wellness Company, Bengaluru	DR ASHOK B. K. DR VIJENDRA PRAKASH ( <i>Alternate</i> )
Indian Medicines Pharmaceutical Corporation Limited, Mohan, Almora	SHRI KAVI RAJ RAI DR BALAJI PANIGRAHI ( <i>Alternate I</i> ) SHRI SRINIWAS CHAUDHARY ( <i>Alternate II</i> )
Indian Institute of Science, Department of Materials Engineering, Bengaluru	PROF RAJEEV RANJAN
Jamia Hamdard, New Delhi	PROF SALEENA KUZHUPPIL BASHIR PROF JAVED ALI ( <i>Alternate I</i> ) DR MD WASI AKHTAR ( <i>Alternate II</i> )
Jamia Millia Islamia, New Delhi	DR JAVID ALI DR NAJMUL ARFIN ( <i>Alternate I</i> ) DR AFREEN INAM ( <i>Alternate II</i> )
National Commission for Indian System of Medicine (NCISM), New Delhi	DR SYED MOHD ABBAS ZAIDI DR AMANULLAH HAJI ( <i>Alternate</i> )
National Institute of Unani Medicine, Bengaluru	PROF ABDUL WADUD DR HAMIDUDDIN ( <i>Alternate I</i> ) DR TASFIYA HAKEEM ANSARI ( <i>Alternate II</i> )
National Medicinal Plant Board, New Delhi	DR R. MURUGESWARAN DR CHINMAY RATH ( <i>Alternate</i> )

<i>Organization</i>	<i>Representative(s)</i>
Pharmacopoeia Commission for Indian Medicine & Homoeopathy, Ghaziabad	SHRI ANUPAM MAURYA SHRI SATISH KUMAR ( <i>Alternate</i> )
Traditional Knowledge Digital Library Unit, Council of Scientific & Industrial Research New Delhi	DR VIJAYLAKSHMI ASTHANA
Unani Drug Manufacturers Association, Delhi	DR NABEEL ANWAR DR SYED MUNEER AZMAT ( <i>Alternate I</i> ) DR WASIF ARBAB ( <i>Alternate II</i> )
BIS Directorate General	SHRIMATI RACHNA SEHGAL, SCIENTIST 'G' AND HEAD (AYUSH) [REPRESENTING DIRECTOR GENERAL ( <i>Ex-officio</i> )]

*Member Secretary*  
DR MAHE ALAM  
SCIENTIST 'C'/DEPUTY DIRECTOR  
(AYUSH), BIS

Panel for Single Herbs Materials and Related Subject, AYD 04/Panel 2

<i>Organization</i>	<i>Representative (s)</i>
Central Council for Research in Unani Medicine, New Delhi	DR R. P. MEENA ( <b><i>Convener</i></b> )
Bureau of India Standards, New Delhi	DR RAGHAVEDERA NAIK
Central Council for Research in Unani Medicine, New Delhi	DR MUSTEHASAN
Central Council for Research in Unani Medicine, New Delhi	DR MOKHTAR ALAM
Central Council for Research in Unani Medicine, New Delhi	SHRIMATI SADAF SUBHANI
Delhi Pharmaceutical Sciences and Research University, New Delhi	PROF AJAY SHARMA
Delhi Pharmaceutical Sciences and Research University, New Delhi	DR ARYA LAKSHMI MARISSETTI
Pharmacopoeia Commission of Indian Medicine and Homoeopathy, (PCIM&H), Ghaziabad	DR ANUPAM MAURYA





## Bureau of Indian Standards

BIS is a statutory institution established under the *Bureau of Indian Standards Act, 2016* to promote harmonious development of the activities of standardization, marking and quality certification of goods and attending to connected matters in the country.

### Copyright

BIS has the copyright of all its publications. No part of these publications may be reproduced in any form without the prior permission in writing of BIS. This does not preclude the free use, in the course of implementing the standard, of necessary details, such as symbols and sizes, type or grade designations. Enquiries relating to copyright be addressed to the Head (Publication & Sales), BIS.

### Review of Indian Standards

Amendments are issued to standards as the need arises on the basis of comments. Standards are also reviewed periodically; a standard along with amendments is reaffirmed when such review indicates that no changes are needed; if the review indicates that changes are needed, it is taken up for revision. Users of Indian Standards should ascertain that they are in possession of the latest amendments or edition by referring to the website-[www.bis.gov.in](http://www.bis.gov.in) or [www.standardsbis.in](http://www.standardsbis.in).

This Indian Standard has been developed from Doc No.: AYD 04 (23340).

### Amendments Issued Since Publication

Amend No.	Date of Issue	Text Affected

## BUREAU OF INDIAN STANDARDS

### Headquarters:

Manak Bhavan, 9 Bahadur Shah Zafar Marg, New Delhi 110002

Telephones: 2323 0131, 2323 3375, 2323 9402

Website: [www.bis.gov.in](http://www.bis.gov.in)

### Regional Offices:

	Telephones
Central : 601/A, Konnectus Tower -1, 6 <sup>th</sup> Floor, DMRC Building, Bhavbhuti Marg, New Delhi 110002	{ 2323 7617
Eastern : 8 <sup>th</sup> Floor, Plot No 7/7 & 7/8, CP Block, Sector V, Salt Lake, Kolkata, West Bengal 700091	{ 2367 0012 2320 9474
Northern : Plot No. 4-A, Sector 27-B, Madhya Marg, Chandigarh 160019	{ 265 9930
Southern : C.I.T. Campus, IV Cross Road, Taramani, Chennai 600113	{ 2254 1442 2254 1216
Western : Manakalya, 4 <sup>th</sup> Floor, NTH Complex (W Sector), F-10, MIDC, Andheri (East), Mumbai 400093	{ 283 25838

**Branches :** AHMEDABAD, BENGALURU, BHOPAL, BHUBANESHWAR, CHANDIGARH, CHENNAI, COIMBATORE, DEHRADUN, DELHI, FARIDABAD, GHAZIABAD, GUWAHATI, HARYANA, HUBLI, HYDERABAD, JAIPUR, JAMMU & KASHMIR, JAMSHEDPUR, KOCHI, KOLKATA, LUCKNOW, MADURAI, MUMBAI, NAGPUR, NOIDA, PARWANOO, PATNA, PUNE, RAIPUR, RAJKOT, SURAT, VIJAYAWADA.