

बिरंज साथी (*Oryza sativa* L.) पारम्परिक औषधियों में उपयोग हेतु सूखे फल- विशिष्टि

Draft Indian Standard

***Biranj Sathī* (*Oryza sativa* L.)- Dried fruit for use in Traditional Medicine –
Specification**

Unani Sectional Committee, AYD-04

Last Date of Comments: 07 December, 2024

FOREWORD

(Formal Clauses will be added later)

This Indian Standard was adopted by Bureau of Indian Standards, after the draft finalized by Unani Sectional Committee had been approved by the Ayush Division Council.

Biranj Sathī consists of dried fruit of *Oryza sativa* L. (Fam. Poaceae); an annual herb, cultivated throughout India. The dried fruits of *Biranj Sathī* are used as an ingredient of different formulations (*Habb-i-Mushil*, *Habb-i-Mulayyin*, *Habb-i Shīr al-Agh*) in Unani medicine.

Biranj Sathī is known by different names in other languages such as Arruz (Arabic); Dhan, Chaul, Chanval (Bengali); Rice, Paddy (English); Bhat, Chorya, Chokha (Gujarati); Chaval, Dhan (Hindi); Akkiege, Nelli (Kannad); Ari (Marathi); Tandul, Sali Bhat (Malayalam); Rambal; Tandulama, Dhanya (Sanskrit); Arshi, Nelli (Tamil); Dhanyamu, Vadlu, Biyyamu (Telugu); Chawal (Urdu).

This standard is one of the series of standards being brought out on ingredients used in formulations of Unani medicine for the advantage of researchers, academicians, students, clinical practitioners and drug manufacturers etc.

In the formulation of this standard, significant assistance has been derived from the Unani Pharmacopoeia of India, Part 1, Vol VII, 2022 published by the Pharmacopoeia Commission for Indian Medicine & Homoeopathy, Ministry of Ayush, Government of India. Inputs have also been derived from the information available in the public domain in print and electronic media including authoritative books.

In the formulation of this standard due consideration has been given to the provisions of *The Drugs and Cosmetics Act*, 1940 and Rules framed thereunder. However, this standard is subject to the restrictions imposed under these Rules and Regulations, wherever applicable.

For the purpose of deciding whether a particular requirement of this standard is complied with, the final value, observed or calculated, expressing the result of a test, shall be rounded off in accordance with **IS 2 : 2022** 'Rules for rounding off numerical values (second revision)'. The number of significant places retained in the rounded off value shall be the same as that of the Specified value in this standard.

BIRANJ SATHĪ (*Oryza sativa* L. -Fruit) DRIED FRUIT FOR USE IN TRADITIONAL MEDICINE — SPECIFICATION

1 SCOPE

This standard prescribes the requirements and methods of test for *Biranj Sathī* which consist of dried fruit of *Oryza sativa* L. (Family Poaceae).

2 REFERENCES

The standards listed in annex-A contain provisions which, through reference in this text, constitute provisions of this standard. At the time of publication, the editions indicated were valid. All standards are subject to revision and parties to agreements based on this standard are encouraged to investigate the possibility of applying the most recent editions of the standards.

3. REQUIREMENTS

3.1 Description:

3.1.0 Macroscopic Examination of *Biranj Sathī* dried fruit

Fruits are small, one seeded, caryopsis, about 0.6-1 cm long and 0.2-0.3 cm wide. They are oblong to ovoid in shape, somewhat angular, blunt and sometimes pointed. The surface is rough due to presence of minute trichomes, faintly longitudinal ridges and furrows which are present in mostly 6 rows and are somewhat compressed, flattened and tightly enclosed by lemma and palea and yellowish-brown in colour. Seed is smooth upto 0.6cm long, oval to oblong and slightly flattened. It is light creamy to white in colour and blunt, oblique, slightly angled in embryo region. Odour is not characteristic and taste is sweet.

3.2.0 Microscopic Examination of *Biranj Sathī* dried fruit

Fruit shows wavy irregular outline. Pericarp and testa are fused together. Pericarp consists of single layered thick lignified sclerenchymatous outer epidermis with clear pits covered by a few thick, blunt, sometimes pointed trichomes and 2-3 layered circular to oval fibre. These are followed by 3-5 layered, tangentially elongated, thick-walled, tabular parenchymatous cells, having a few scattered fibro vascular bundles and single layered, thin, elongated, slightly wavy inner epidermal cells. Testa consists of thin walled, elongated, 2-3 layered parenchymatous cells with interrupted tube cells followed by single layered, oval to rectangular, parenchymatous layer containing aleurone grains. Endosperm is albuminous, consisting of wide, thin-walled, elongated to polygonal, parenchymatous cells packed with numerous, minute, single polyhedral starch grains, having hilum without concentric striations, measuring 3-12 μ in diameter and compound starch grains 2-150 components. Embryo is small, lying in a groove at one end of the endosperm and separated by a layer of epithelium. Embryo consists of a shield shaped cotyledon known as scutellum.

3.3 Macroscopic and Microscopic Examination of Powder of *Biranj Sathī* dried fruit

Powder is light cream in colour. Odour is not characteristic. Taste is slightly sweet. The microscopy shows fragments of elongated thick-walled, lignified sclerenchymatous epidermal cells with clear pits. Trichomes are unicellular with pointed tips upto 700 μ . Parenchyma cells are oval to rectangular thin-walled cells filled with aleurone grains. Endosperm contains numerous, thin-walled, elongated to polygonal, parenchyma cells packed with numerous, minute, single polyhedral starch grains. Starch grains are simple measuring upto 12 μ and compound starch grains (2-150 components) (Fig.1).

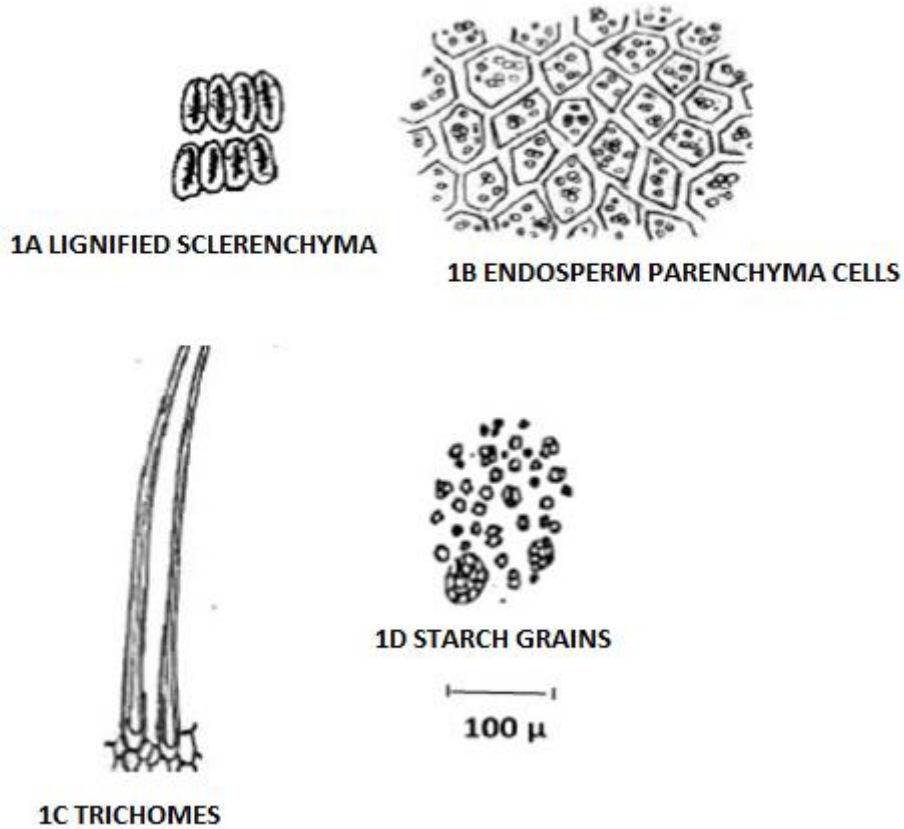


Fig. 1. Biranj Sathī (*Oryza sativa* L.-Fruit) - Powder Microscopy

4.1 General

4.1.1 *Biranj Sathī* (*Oryza sativa* L.-Fruit) shall be free from extraneous/ artificial flavours.

4.1.2 *Biranj Sathī* (*Oryza sativa* L.-Fruit) shall comply with physical, chemical and microbiological requirements given in Table 1.

4.1.3 If any pesticides other than those for which minimum requirements are given in Table 1 are applied to the herb before or after harvesting, those should also be tested. Their limit shall be calculated using the following formula:

$$\frac{ADI \times M}{MDD \times 100}$$

where

ADI = Acceptable daily intake of pesticide as published by Food and Agriculture Organization-World Health Organization (FAO-WHO), in milligrams per kilogram of body mass,

M = Body mass in kilograms (60 kg).

MDD = Maximum daily dose of the drug, in kilograms.

5 PACKING, STORAGE AND MARKING

5.1 Packing

Biranj Sathī (*Oryza sativa* L.-Fruit) shall be packed in clean, sound and dry container made up of glass, foodgrade polymers, wood or jute bags. The wooden boxes or jute bags shall be suitably lined with moisture proof lining which does not impart any foreign smell to the product. The packing material shall be free from any fungal or insect infestation and should not impart any foreign smell. Each container shall be securely closed and sealed.

5.2 Storage

Biranj Sathī (*Oryza sativa* L.-Fruit) shall be stored under conditions that prevent contamination and, as far as possible, deterioration. Storage area shall be clean, well ventilated, protected from light, moisture, insects and rodents.

5.3 Marking

The following particulars shall be legibly and indelibly marked or labelled on each pack:

- a) Name of the material including part of the plant, botanical name, and trade name or brand name, if any;
- b) Name and address of the producer or packer;
- c) Batch number;
- d) Net quantity;
- e) State and country of production;
- f) Date of packing (MM/YYYY);
- g) Instructions for storage;
- h) Any other information requested by the buyer, such as the date of harvesting (MM/YYYY) (if known).

The above information, or part of it, may instead appear in the documentation after agreement between the buyer and the seller.

5.4 BIS Certification Marking

The product(s) conforming to the requirements of this standard may be certified as per the conformity assessment schemes under the provisions of the *Bureau of Indian Standards Act, 2016* and the Rules and Regulations framed there under, and the product(s) may be marked with the Standard Mark.

6 SAMPLING

6.1 Representative samples of the material shall be drawn and tested for conformity to this specification as prescribed in IS 13145.

6.2 The samples of *Biranj Sathī* dried fruit shall be tested to ascertain conformity of the material to the requirements in accordance with the relevant clauses given in column 4 of Table 1.

WC DRAFT DOCUMENT: FOR COMMENTS ONLY
(Not to be reproduced without permission of BIS or used as Indian Standard)

AYD-04 (26483)
October, 2024

Table 1 Requirements for *Biranj Sathī* (*Oryza sativa* L.-Fruit)
(Clauses 3.2.2, 3.2.3 and 5.2)

Sl No.	Characteristic	Requirement	Method of Test, Ref to
(1)	(2)	(3)	(4)
i)	Foreign matter, percent by mass, <i>Max</i>	<2%	6.2 of IS 4333 (Part 1)
ii)	Total ash, percent by mass, <i>Max</i>	< 6%	6 of IS 1797
iii)	Acid insoluble ash, percent by mass, <i>Max</i>	< 5 %	8 of IS 1797
iv)	Alcohol soluble Extractive, percent by mass, <i>Min</i>	> 1%	10 of IS 1797
v)	Water soluble Extractive, percent by mass, <i>Min</i>	> 1 %	11 of IS 1797
vi)	Thin-layer chromatogram identification	-	Annex B
vii)	Aflatoxin (B1), ng/kg, <i>Max</i>	2.0	IS 16287
viii)	Aflatoxin (B1+B2+G1+G2), ng/kg, <i>Max</i>	5.0	IS 16287
ix)	Lead, mg/kg, <i>Max</i>	10.0	IS 16913
x)	Mercury, mg/kg, <i>Max</i>	1.0	IS 16913
xi)	Cadmium, mg/kg, <i>Max</i>	0.3	IS 16913
xii)	Arsenic, mg/kg, <i>Max</i>	3.0	IS 16913
xiii)	<i>Staphylococcus aureus</i> , per g	Absent for extract and powder	Annex C
		Absent for plant material	
xiv)	<i>Pseudomonas aeruginosa</i> , per g	Absent for extract and powder	Annex C
		Absent for plant material	
xv)	<i>Salmonella</i> spp., per g	Absent	Annex C
xvi)	<i>Escherichia coli</i> , per g	Absent for extract and powder	Annex C
		10 for plant material	
xvii)	Total microbial plate count, per g, <i>Max</i>	10 ⁵ for extract and powder	Annex C
		10 ⁷ for plant material	
xviii)	Total yeast and mould, per g, <i>Max</i>	10 ³ for extract and powder	Annex C

WC DRAFT DOCUMENT: FOR COMMENTS ONLY
(Not to be reproduced without permission of BIS or used as Indian Standard)

AYD-04 (26483) October, 2024

		10 ⁵ for plant material	
xix)	Pesticide residues		
1.	Alachlor, mg/kg, <i>Max</i>	0.02	
2.	Aldrin and dieldrin (sum of), mg/kg, <i>Max</i>	0.05	
3.	Azinphos-methyl, mg/kg, <i>Max</i>	1.0	
4.	Bromopropylate, mg/kg, <i>Max</i>	3.0	
5.	Chlordane (sum of cis-, trans - and Oxythlordane), mg/kg, <i>Max</i>	0.05	
6.	Chlorfenvinphos, mg/kg, <i>Max</i>	0.5	
7.	Chlorpyrifos, mg/kg, <i>Max</i>	0.2	
8.	Chlorpyrifos-methyl, mg/kg, <i>Max</i>	0.1	
9.	Cypermethrin (and isomers), mg/kg, <i>Max</i>	1.0	
10.	DDT (sum of p,p'-DDT, o,p'-DDT, p,p'-DDE and p,p'-TDE), mg/kg, <i>Max</i>	1.0	IS 17924
11.	Deltamethrin, mg/kg, <i>Max</i>	0.5	
12.	Diazinon, mg/kg, <i>Max</i>	0.5	
13.	Dichlorvos, mg/kg, <i>Max</i>	1.0	
14.	Dithiocarbamates (as CS ₂), mg/kg, <i>Max</i>	2.0	
15.	Endosulfan (sum of isomers and endosulfan sulphate), mg/kg, <i>Max</i>	3.0	
16.	Endrin, mg/kg, <i>Max</i>	0.05	
17.	Ethion, mg/kg, <i>Max</i>	2.0	
18.	Fenitrothion, mg/kg, <i>Max</i>	0.5	
19.	Fenvalerate, mg/kg, <i>Max</i>	1.5	
20.	Fonofos, <i>Max</i>	0.05	
21.	Heptachlor (sum of heptachlor and heptachlor epoxide), mg/kg, <i>Max</i>	0.05	
22.	Hexachlorobenzene, mg/kg, <i>Max</i>	0.1	
23.	Hexachlorocyclohexane isomers (other than γ), mg/kg, <i>Max</i>	0.3	
24.	Lindane (γ -hexachlorocyclohexane), mg/kg, <i>Max</i>	0.6	
25.	Malathion, mg/kg, <i>Max</i>	1.0	
26.	Methidathion, mg/kg, <i>Max</i>	0.2	
27.	Parathion, mg/kg, <i>Max</i>	0.5	IS 17924
28.	Parathion-methyl, mg/kg, <i>Max</i>	0.2	
29.	Permethrin, mg/kg, <i>Max</i>	1.0	
30.	Phosalone, mg/kg, <i>Max</i>	0.1	
31.	Piperonyl butoxide, mg/kg, <i>Max</i>	3.0	
32.	Pirimiphos-methyl, mg/kg, <i>Max</i>	4.0	
33.	Pyrethrins (sum of), mg/kg, <i>Max</i>	3.0	

WC DRAFT DOCUMENT: FOR COMMENTS ONLY
(Not to be reproduced without permission of BIS or used as Indian Standard)

AYD-04 (26483)
October, 2024

34.	Quintozene (sum of quintozene, pentachloroaniline and methyl pentachlorophenyl sulphide), mg/kg, <i>Max</i>	1.0	
-----	---	-----	--

7 QUALITY OF REAGENTS

7.1 Reagents including pure chemicals used shall be of analytical grade.

7.2 Reagent grade water for laboratory use shall be as per IS 1070.

NOTE – ‘Pure chemicals’ shall mean chemicals that do not contain impurities which effect the results of analysis.

ANNEX-A

LIST OF STANDARDS REFERRED

(Clause 2)

<i>IS No.</i>	<i>Title</i>
1070 : 2023	Reagent grade water – Specification (<i>fourth revision</i>)
1797 : 2017	Spices and condiments – Methods of test (<i>third revision</i>)
4333 : 2018	Methods of analysis for foodgrains: Part 1 Refractions (<i>third revision</i>)
13145 : 2014	Spices and condiments – Methods of sampling (<i>second revision</i>)
16913 : 2018	Methods of test for cosmetics – Determination of heavy metals (Arsenic, Cadmium, Lead and Mercury) by atomic absorption spectrometry (AAS)
17924 : 2022	Determination of pesticide residue in herbal materials
16287: 2015/ ISO 16050 : 2003	Foodstuffs – Determination of aflatoxin B ₁ , and the total content of aflatoxins B ₁ , B ₂ , G ₁ and G ₂ in cereals, nuts and derived products – High performance liquid chromatographic method

ANNEX-B

[Table 1, SI No. vi), Col (4)]

Thin-layer chromatogram identification

GENERAL

The chromatographic profile of the test solution is similar to that of the reference solution, prepared with the BRS *Biranj Sathi* (*Oryza sativa* L.-Fruit) under the same conditions, with respect to the position and fluorescence/colour of the bands.

B-1 THIN LAYER CHROMATOGRAPHY (TLC)

B-2.1 Apparatus

B-2.1.1 *Thin Layer Chromatograph with a Pre Coated Silica Gel 60F₂₅₄ Plate*

B-2.1.2 *Analytical Balance*

B-2.1.3 *Standard Glassware*

B-2.1.4 *Water Bath*

B-2.1.5 *TLC plate development Chamber with twin truff*

B-2.1.5 *Hot air oven*

B-2.1.5 *Spraying apparatus*

B-2.2 Reagents

B-2.2.1 *Preparation of test Solution*

B-2.2.1 *Reference Standard – Botanical Reference Standard*

B-2.2.2 *Chloroform – AR or Equivalent Grade*

B-2.2.3 *Ethanol – AR or Equivalent Grade*

B-2.2.4 *Toluene – AR or Equivalent Grade*

B-2.2.5 *Ethyl acetate – AR or Equivalent Grade*

B-2.2.6 *Formic acid – AR or Equivalent Grade*

B-2.2.7 *Vanillin-sulphuric acid reagent – AR or Equivalent Grade*

B-2.3 Procedure

B-2.3.1 *Solvent System* Dissolve chloroform, Alcohol and toluene, ethyl acetate and formic acid in the ratio (8.7:1.3:0.1)

B-2.3.2 *Preparation of Test Solution* To 2 g of the substance being examined, add 20 ml of methanol, heat on a water bath for 10 to 15 min, cool and filter.

B-2.2.1 *Preparation of Botanical Reference Standard (BRS) Solution.*

B 3.0 Identification:

B 3.1 *Thin Layer Chromatography of Chloroform Extract:*

B-3.2 *Thin-layer chromatograph with a Pre-Coated Silica Gel 60F₂₅₄ plate.*

B-3.3 *Test Solution:* Prepare the extract from 2 g of coarsely powdered plant material by refluxing with 20 ml of *chloroform* for 30 min. Filter, concentrate the extract and use as test solution.

B-3.4 **Procedure:** Apply 10 µl of the test solution as 8 mm band at a height of 10 mm from the base of a 10 x 10 cm TLC plate. Develop the plate to a distance of 8 cm from the band using solvent system: Toluene: Ethyl acetate: Formic acid (8.7:1.3:0.1). Air dry the plate and examine under UV 254 nm and UV 366 nm. Spray the plate with Vanillin-sulphuric acid reagent and heat at temperature of 105-110°C till the bands appear.

The chromatographic profile of the test solution is similar to that of the reference solution, prepared with the BRS *Biranj Sathi* (fruit) under the same conditions, with respect to the position and fluorescence/colour of the bands (Fig.2).

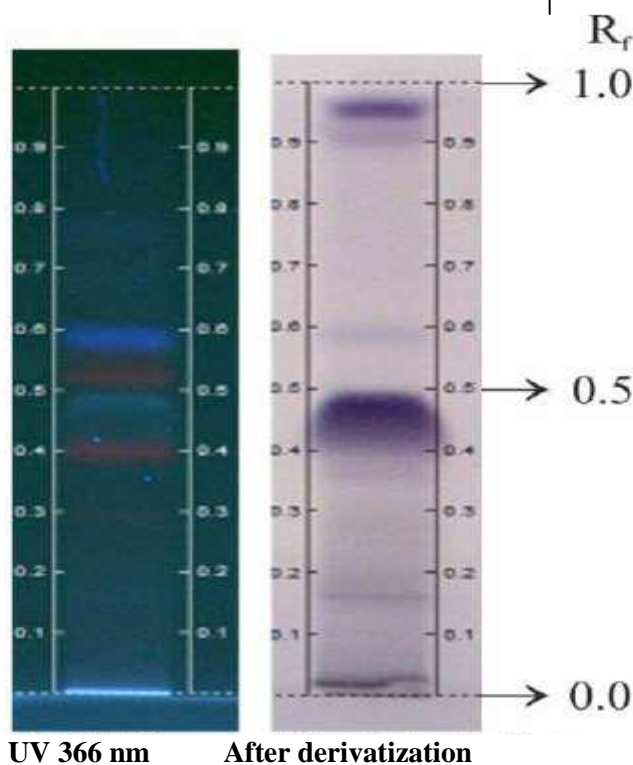


Fig. 2: TLC fingerprints of *Biranj Sathi* (*Oryza sativa* L.-Fruit)

B-4.0 *Thin Layer Chromatography of Alcohol Extract:*


B-4.1 *Thin-layer chromatograph with a Pre-Coated Silica Gel 60F₂₅₄ plate.*

B-4.2 *Test Solution:* Prepare the extract from 2 g of coarsely powdered plant material by refluxing with 20ml of *ethanol* for 30 min. Filter, concentrate the extract and use as test solution.

B-4.3 *Procedure:* Apply 10 µl of the test solution as 8 mm band at a height of 10 mm from the base of a 10 x 10 cm TLC plate. Develop the plate to a distance of 8 cm from the band using solvent system: *toluene: ethyl acetate: formic acid* (8:2:0.1). Air dry the plate and examine under UV 366 nm. Spray the plate with *vanillin-sulphuric acid reagent* and heat at temperature of 105-110°C till the bands appear.

The chromatographic profile of the test solution is similar to that of the reference solution, prepared with the BRS *Biranj Sathi* (fruit) under the same conditions, with respect to the position and fluorescence/colour of the bands (Fig. 3).

AYD-04 (26483)
October, 2024

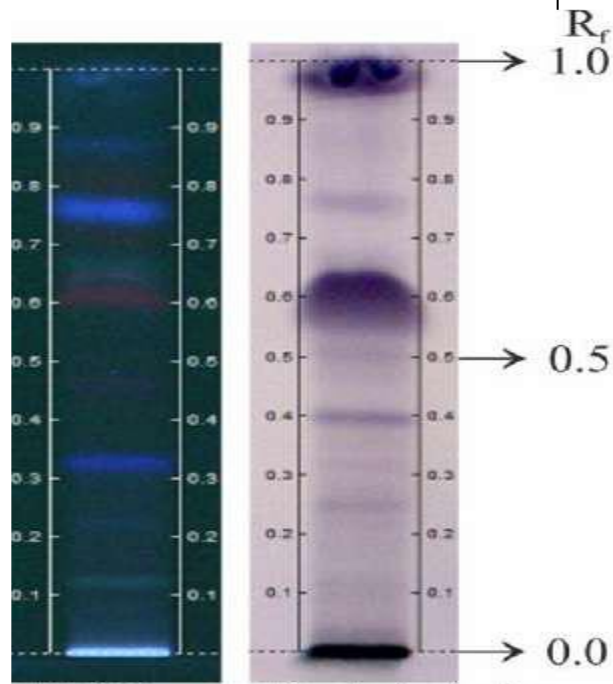

UV 366 nm After derivatization

Fig. 3: TLC fingerprints of *Biranj Sathī* (*Oryza sativa* L.-Fruit)

ANNEX C
TEST FOR MICROBIAL LIMITS

[Table 1, Sl No. (xii), (xiii), (xiv), (xv), (xvi), and (xvii)]

C-1 GENERAL

The tests are designed for the estimation of the number of viable aerobic microorganisms present and for detecting the presence of designated microbial species in the extract. The term 'growth' is used to designate the presence and presumed proliferation of viable microorganisms.

C-2 APPARATUS

C-2.1 Oven for Dry Sterilization

C-2.2 Autoclave for Wet Sterilization

C-2.3 Incubator

C-2.4 Water Bath

C-2.5 pH-meter

C-2.6 Sterile Membrane Filters, 50 mm in Diameter

C-2.7 Colony Counting Equipment

C-2.8 Analytical Balance

C-2.9 Standard Glassware

C-3 MEDIA

Culture media may be prepared as given below or dehydrated culture media may be used provided that, when reconstituted as directed by the manufacturer, they have similar ingredients and / or yield media comparable to those obtained from the formulae given below.

Where agar is specified in a formula, use agar that has a moisture content of not more than 15 percent. Where water is called for in a formula, use purified water. Unless otherwise indicated, the media should be sterilized by heating in an autoclave (15 psi) at 121 °C for 15 min. In preparing media by the formulas given below, dissolve the soluble solids in the water, using heat, if necessary, to effect complete solution, add solutions of 0.1N hydrochloric acid or 0.1N sodium hydroxide in quantities sufficient to yield the required pH in the medium when it is ready for use. Determine the pH at 25 °C ± 2 °C.

C-3.1 Baird Parker Agar Medium

C-3.1.1 Composition

Pancreatic digest of casein	10.0 g
Beef extract	5.0 g
Yeast extract	1.0 g
Lithium chloride	5.0 g
Agar	20.0 g
Glycine	12.0 g
Sodium pyruvate	10.0 g
Water	1 000 ml

C-3.1.2 Preparation

Suspend the components in 1 000 ml of water, heat with frequent agitation and boil for 1 min. Sterilize, cool in between 45 °C to 50 °C, add 10 ml of a one percent w/v solution of sterile potassium tellurite and 50 ml of egg yolk emulsion. Mix thoroughly, but gently and pour into plates. (Prepare the egg yolk emulsion by disinfecting the surface of whole shell eggs, aseptically cracking the eggs, and separating out intact yolks into a sterile graduated cylinder. Add sterile saline solution, get a 3 to 7 ratio of egg yolk to saline. Add to a sterile blender cup and mix at high speed for 5 seconds). Adjust the pH after sterilization to 6.8 ± 0.2 .

C-3.2 Bismuth Sulphite Agar Medium

C-3.2.1 Composition

Solution (1)

Beef extract	6 g
Peptone	10 g
Agar	24 g
Ferric citrate	0.4 g
Brilliant green	10 mg
Water	1 000 ml

Solution (2)

Ammonium bismuth citrate	3 g
Sodium sulphite	10 g
Anhydrous disodium hydrogen phosphate	5 g
Dextrose monohydrate	5 g
Water	100 ml

C-3.2.2 Preparation

Suspend the components of Solution 1 in 1 000 ml of water. Heat to boiling to dissolve the medium completely. Sterilize by maintaining at 115 °C for 30 min.

Suspend the components of Solution 2 in 100 ml of water. Heat to boiling to dissolve the medium completely. Do not autoclave.

Add 1 volume of solution 2 to 10 volumes of solution 1 previously melted and cooled to a temperature of 55 °C. Bismuth Sulphite Agar Medium should be stored at 2 °C to 8 °C for 5 days before use.

C-3.3 Brilliant Green Agar Medium

C-3.3.1 Composition

Peptone	10.0 g
Yeast extract	3.0 g
Lactose	10.0 g
Sucrose	10.0 g
Sodium chloride	5.0 g
Phenol red	80.0 g
Brilliant green	12.5 mg
Agar	12.0 g
Water	1 000 ml

C-3.3.2 Preparation

Mix the components, allow to stand for 15 min, sterilize by maintaining at 115 °C for 30 min and mix before pouring.

C-3.4 Buffered Sodium Chloride Peptone Solution pH 7.0

C-3.4.1 Composition

Potassium dihydrogen phosphate	3.56 g
Disodium hydrogen phosphate	7.23 g
Sodium chloride	4.30 g
Peptone (meat or casein)	1.0 g
Water	1 000 ml

C-3.4.2 Preparation

Mix the components and heat if necessary to dissolve the medium completely. 0.1 to 1.0 percent w/v Polysorbate 20 or polysorbate 80 may be added. Sterilize by heating in an autoclave at 121 °C for 15 min.

C-3.5 Casein Soybean Digest Agar Medium

C-3.5.1 Composition

Pancreatic digest of casein	15.0 g
Papaic digest of soybean meal	5.0 g
Sodium chloride	5.0 g
Agar	15.0 g
Water	1 000 ml

C-3.5.2 Preparation

Suspend the components in water. Heat to boiling to dissolve the medium completely. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to 7.3 ± 0.2 .

C-3.6 Cetrinide Agar Medium

C-3.6.1 Composition

Pancreatic digest of gelatin	20.0 g
Magnesium chloride	1.4 g
Potassium sulphate	10.0 g
Cetrinide	0.3 g
Agar	13.6 g
Glycerin	10.0 g
Water	1 000 ml

C-3.6.2 Preparation

Dissolve the components in 1 000 ml of water. Heat to boiling for 1 min with shaking. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to 7.0 to 7.4.

C-3.7 Deoxycholate Citrate Agar Medium

C-3.7.1 Composition

Beef extract	5.0 g
Peptone	5.0 g
Lactose	10.0 g
Trisodium citrate	8.5 g
Sodium thiosulphate	5.4 g
Ferric citrate	1.0 g
Sodium deoxycholate	5.0 g
Neutral red	0.02 g
Agar	12.0 g
Water	1 000 ml

C-3.7.2 Preparation

Mix the components and allow to stand for 15 min. Gently boil with continuous stirring and continue boiling until solution is complete. Cool to 80 °C, mix, pour and cool rapidly.

Care should be taken not to overheat Deoxycholate Citrate Agar during preparation. It should not be re-melted and the surface of the plates should be dried before use.

C-3.8 Fluid Casein Digest Soya Lecithin Polysorbate 20 Medium

C-3.8.1 Composition

Pancreatic digest of casein	20 g
Soya lecithin	5 g
Polysorbate 20	40 ml
Water	1 000 ml

C-3.8.2 Preparation

Dissolve the pancreatic digest of casein and soya lecithin in water, heating in a water bath at 48 °C to 50 °C for about 30 min to effect solution. Add polysorbate 20, mix and dispense as desired. Sterilize at 121 °C for 15 min in an autoclave.

C-3.9 Fluid Lactose Medium

C-3.9.1 Composition

Beef extract	3.0 g
Pancreatic digest of gelatin	5.0 g
Lactose	5.0 g
Water	1 000 ml

C-3.9.2 Preparation

Suspend the components in 1 000 ml water. Heat if necessary to dissolve the medium completely. Sterilize at 121 °C for 15 min in an autoclave. Cool as quickly as possible after sterilization. Adjust the pH after sterilization to 6.9 ± 0.2.

C-3.10 Lactose Broth Medium

C-3.10.1 Composition

Beef extract	3.0 g
Pancreatic digest of gelatin	5.0 g
Lactose	5.0 g
Water	1 000 ml

C-3.10.2 Preparation

Suspend the components in water and heat if necessary to dissolve the medium completely. Sterilize at 121 °C for 15 min in an autoclave. Adjust the pH after sterilization to 6.9 ± 0.2 .

C-3.11 Levine Eosin Methylene Blue Agar Medium

C-3.11.1 Composition

Pancreatic digest of gelatin	10.0 g
Dibasic potassium phosphate	2.0 g
Agar	15.0 g
Lactose	10.0 g
Eosin Y	400 mg
Methylene blue	65 mg
Water	1 000 ml

C-3.11.2 Preparation

Dissolve the pancreatic digest of gelatin, dibasic potassium phosphate and agar in water with warming and allow to cool. Just prior to use, liquify the gelled agar solution and the remaining ingredients, as solutions, in the following amounts and mix. For each 100 ml of the liquified agar solution use 5 ml of a 20 percent w/v solution of lactose, 2 ml of a 2 percent w/v solution of eosin Y and 2 ml of a 0.33 percent w/v solution of methylene blue. The finished medium may not be clear. Adjust the pH after sterilization to 7.1 ± 0.2 .

C-3.12 MacConkey Agar Medium

C-3.12.1 Composition

Pancreatic digest of gelatin	17.0 g
Peptone (meat and casein, equal parts)	3.0 g
Lactose	10.0 g
Sodium chloride	5.0 g
Bile salts	1.5 g
Agar	13.5 g
Neutral red	30 mg
Crystal violet	1 mg
Water	1 000 ml

C-3.12.2 Preparation

Boil the mixture of solids and water for 1 min to effect solution. Sterilize at 121 °C for 15 min in an autoclave. Adjust the pH after sterilization to 7.1 ± 0.2 .

C-3.13 MacConkey Broth Medium

C-3.13.1 Composition

Pancreatic digest of gelatin	20.0 g
Lactose	10.0 g
Dehydrated ox bile	5.0 g
Bromocresol purple	10 mg
Water	1 000 ml

C-3.13.2 Preparation

Suspend the components in 1 000 ml of water and heat if necessary to dissolve the medium completely. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to 7.3 ± 0.2 .

C-3.14 Mannitol Salt Agar Medium

C-3.14.1 Composition

Pancreatic digest of gelatin	5.0 g
Peptic digest of animal tissue	5.0 g
Beef extract	1.0 g
D-Mannitol	10.0 g
Sodium chloride	75.0 g
Agar	15.0 g
Phenol red	25 mg
Water	1 000 ml

C-3.14.2 Preparation

Mix the components, heat with frequent agitation and boil for 1 min to effect solution. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to 7.4 ± 0.2 .

C-3.15 Nutrient Broth Medium

C-3.15.1 Composition

Beef extract	10.0 g
Peptone	10.0 g
Sodium chloride	5 mg
Water	1 000 ml

C-3.15.2 Preparation

Dissolve the components with the aid of heat. Adjust the pH to 8.0 to 8.4 with 5M sodium hydroxide and boil for 10 min. Filter and sterilize by maintaining at 115 °C for 30 min and adjust pH to 7.3 ± 0.1 .

C-3.16 Nutrient Agar Medium

Nutrient broth gelled by the addition of 1 to 2 percent w/v of agar.

C-3.17 Pseudomonas Agar Medium for Detection of Fluorescein

C-3.17.1 Composition

Pancreatic digest of casein	10.0 g
Peptic digest of animal tissue	10.0 g
Anhydrous dibasic potassium phosphate	1.5 g
Magnesium sulphate hepta hydrate	1.5 g
Glycerin	10.0 ml
Agar	15.0 g
Water	1 000 ml

C-3.17.2 Preparation

Dissolve the solid components in water before adding glycerin. Heat with frequent agitation and boil for 1 min to effect solution. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to 7.2 ± 0.2 .

C-3.18 Pseudomonas Agar Medium for Detection of Pyocyanin

C-3.18.1 Composition

Pancreatic digest of gelatin	20.0 g
Anhydrous magnesium chloride	1.4 g
Anhydrous potassium sulphate	10.0 g
Agar	15.0 g
Glycerin	10.0 ml
Water	1 000 ml

C-3.18.2 Preparation

Dissolve the solid components in water before adding glycerin. Heat with frequent agitation and boil for 1 min to effect solution. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to 7.2 ± 0.2 .

C-3.19 Sabouraud Dextrose Agar Medium

C-3.19.1 Composition

Dextrose	40 g
Peptic digest of animal tissue and Pancreatic digest of casein (1:1)	10 g
Agar	15 g
Water	1 000 ml

C-3.19.2 Preparation

Mix the components and heat to boiling to dissolve completely. Sterilize at 121 °C for 15 min in an autoclave and adjust the pH after sterilization to 5.6 ± 0.2 .

C-3.20 Sabouraud Dextrose Agar Medium with Antibiotics

To 1 000 ml of Sabouraud Dextrose Agar Medium, add 0.1 g of benzylpenicillin sodium and 0.1 g of tetracycline HCL or alternatively add 50 mg of chloramphenicol immediately before use.

C-3.21 Selenite F Broth

C-3.21.1 Composition

Peptone	5 g
Lactose	4 g
Disodium hydrogen phosphate	10 g
Sodium hydrogen selenite	4 g
Water	1 000 ml

C-3.21.2 Preparation

Suspend the components in water and mix well. Warm to dissolve the medium completely. Distribute in sterile containers and sterilize by maintaining at 100 °C for 30 min.

C-3.22 Fluid Selenite Cystine Medium

C-3.22.1 Composition

Pancreatic digest of casein	5.0 g
Lactose	4.0 g
Sodium phosphate	10.0 g
Sodium hydrogen selenite	4.0 g
L-Cystine	10.0 mg
Water	1 000 ml

C-3.22.2 Preparation

Suspend the components in water and heat in flowing steam for 15 min. Adjust the final pH to 7.0 ± 0.2. Do not sterilize.

C-3.23 Tetrathionate Broth Medium

C-3.23.1 Composition

Beef extract	0.9 g
Peptone	4.5 g
Yeast extract	1.8 g
Sodium chloride	4.5 g
Calcium carbonate	25.0 g
Sodium thiosulphate	40.7 g
Water	1 000 ml

C-3.23.2 Preparation

Dissolve the solids in 1 000 ml of water and heat the solution to boil. On the day of use, add a solution prepared by dissolving 5 g of potassium iodide and 6 g of iodine in 20 ml of water.

C-3.24 Tetrathionate Bile Brilliant Green Broth Medium

C-3.24.1 Composition

Peptone	8.6 g
Dehydrated ox bile	8.0 g
Sodium chloride	6.4 g
Calcium carbonate	20.0 g
Potassium tetrathionate	20.0 g
Brilliant green	70 mg
Water	1 000 ml

C-3.24.2 Preparation

Suspend the components in 1 000 ml of water. Heat just to boiling. Do not autoclave or reheat. Adjust the pH so that after heating it is 7.0 ± 0.2 .

C-3.25 Triple Sugar Iron Agar Medium

C-3.25.1 Composition

Beef extract	3.0 g
Yeast extract	3.0 g
Peptone	20.0 g
Lactose	10.0 g
Sucrose	10.0 g
Dextrose monohydrate	1.0 g
Ferrous sulphate	0.2 g
Sodium chloride	5.0 g
Sodium thiosulphate	0.3 g
Phenol red	24 mg
Water	1 000 ml

C-3.25.2 Preparation

Mix the components in 1 000 ml of water, allow standing for 15 min. Heat to boiling to dissolve the medium completely. Mix well and distribute into test tubes and Sterilize by maintaining at $121\text{ }^{\circ}\text{C}$ for 15 min. Allow the medium to stand in a sloped form with a butt about 2.5 cm long.

C-3.26 Urea Broth Medium

C-3.26.1 Composition

Potassium dihydrogen orthophosphate	9.1 g
Anhydrous disodium hydrogen phosphate	9.5 g
Urea	20.0 g
Yeast extract	0.1 g
Phenol red	10 mg
Water	1 000 ml

C-3.26.2 Preparation

Mix the components, sterilize by filtration and distribute aseptically in sterile containers.

C-3.27 Vogel Johnson Agar Medium

C-3.27.1 Composition

WC DRAFT DOCUMENT: FOR COMMENTS ONLY
(Not to be reproduced without permission of BIS or used as Indian Standard)

AYD-04 (26483)
October, 2024

Pancreatic digest of casein	10.0 g
Yeast extract	5.0 g
Mannitol	10.0 g
Dibasic potassium phosphate	5.0 g
Lithium chloride	5.0 g
Glycerin	10.0 g
Agar	16.0 g
Phenol red	25.0 mg
Water	1 000 ml

C-3.27.2 Preparation

Suspend the components in 1 000 ml of water. Boil the solution of solids for 1 min. Sterilize at 121 °C for 15 min in an autoclave. Cool to 45 °C to 50 °C and add 20 ml of 1 percent w/v sterile solution of potassium tellurite. Adjust the pH after sterilization to 7.0 ± 0.2 .

C-3.28 Xylose Lysine Deoxycholate Agar Medium

C-3.28.1 Composition

Xylose	3.5 g
l-Lysine	5.0 g
Lactose	7.5 g
Sucrose	7.5 g
Sodium chloride	5.0 g
Yeast extract	3.0 g
Phenol red	80 mg
Agar	13.5 g
Sodium deoxycholate	2.5 g
Sodium thiosulphate	6.8 g
Ferric ammonium citrate	800 mg
Water	1 000 ml

C-3.28.2 Preparation

Suspend the components in 1 000 ml of water. Heat with frequent agitation until the medium boils. Do not overheat or sterilize. Transfer at once to a water bath maintained at about 50 °C and pour into plates as soon as the medium has cooled. Adjust the final pH to 7.4 ± 0.2 .

C-4 SAMPLING

Use 10 ml or 10 g specimens for each of the tests specified in the Indian Standard.

C-5 PRECAUTIONS

The microbial limit tests should be carried out under conditions designed to avoid accidental contamination during the test. The precautions taken to avoid contamination must be such that they do not adversely affect any microorganisms that should be revealed in the test.

C-6 PRELIMINARY TESTING

Prior to doing the tests, inoculate diluted extracts being examined with separate viable cultures of *Escherichia coli*, *Salmonella* species, *Pseudomonas aeruginosa* and *Staphylococcus aureus*. This is done by adding 1 ml of 24 h broth culture containing not less than 1 000 microorganisms to the first

dilution (in buffer solution pH 7.2, fluid soybean casein digest medium or fluid lactose medium) of the test material and following the test procedure. If the organisms fail to grow in the relevant medium the procedure should be modified by (a) increasing the volume of diluent with the quantity of test material remaining the same, or (b) incorporating a sufficient quantity of a suitable inactivating agent in the diluents, or (c) combining the afore mentioned modifications so as to permit growth of the organisms in the media. If inhibitory substances are present in the extracts, 0.5 percent of soya lecithin and 4 percent of polysorbate 20 may be added to the culture medium. Alternatively, repeat the test as described in the previous paragraph, using fluid casein digest soya lecithin polysorbate 20 medium, to demonstrate neutralization of preservatives or other antimicrobial agents in the test material. Where inhibitory substances are contained in the extracts and the latter is soluble, the membrane filtration method described under Total Aerobic Microbial Count may be used. If in spite of incorporation of suitable inactivating agents and a substantial increase in the volume of diluent it is still not possible to recover the viable cultures described above and where the article is not suitable for applying the membrane filtration method, it can be assumed that the failure to isolate the inoculated organism may be due to the bactericidal activity of the product. This may indicate that the article is not likely to be contaminated with the given species of microorganisms. However, monitoring should be continued to establish the spectrum of inhibition and bactericidal activity of the article.

C-7 TOTAL AEROBIC MICROBIAL COUNT

Pre-treat the extracts and raw materials being examined as described below.

NOTE – The raw materials need to be ground as a coarse powder before analysis.

C-7.1 Water Soluble Products

Dissolve 10 g or dilute 10 ml of the extract preparation being examined, unless otherwise specified, in buffered sodium chloride peptone solution pH 7.0 or any other suitable medium shown to have no antimicrobial activity under the conditions of test and adjust the volume to 100 ml with the same medium. If necessary, adjust the pH to about 7.

C-7.2 Products Insoluble in Water (Non Fatty)

Suspend 10 g or 10 ml of the extract preparation being examined, unless otherwise specified, in buffered sodium chloride peptone solution pH 7.0 or any other suitable medium shown not to have antimicrobial activity under the conditions of the test and dilute to 100 ml with the same medium. If necessary, divide the preparation being examined and homogenize the suspension mechanically. A suitable surface active agent such as 0.1 percent w/v of polysorbate 80 may be added to assist the suspension of poorly wettable substances. If necessary, adjust pH of the suspension to about 7.

C-7.3 Fatty Products

Homogenise 10 g or 10 ml of the extract preparation being examined, unless otherwise specified, with 5 g of polysorbate 20 or polysorbate 80. If necessary, heat to not more than 40 °C. Mix carefully while maintaining the temperature in the water bath or in an oven. Add 85 ml of buffered sodium chloride peptone solution pH 7.0 or any other suitable medium shown to have no antimicrobial activity under the conditions of the test, heated to not more than 40 °C if necessary. Maintain this temperature for the shortest time necessary for formation of an emulsion and in any case for not more than 30 min. If necessary, adjust the pH to about 7.

C-7.4 Examination of the Sample

Determine the total aerobic microbial count in the extract being examined by any of the following methods.

C-7.5 Membrane Filtration

Use membrane filters 50 mm in diameter and having a nominal pore size not greater than 0.45 μm the effectiveness of which in retaining bacteria has been established for the type of preparation being examined. Transfer 10 ml or a quantity of each dilution containing 1g of the preparation being examined to each of two membrane filters and filter immediately. If necessary, dilute the pretreated extract preparation so that a colony count of 10 to 100 may be expected. Wash each membrane by filtering through it three or more successive quantities, each of about 100 ml, of a suitable liquid such as buffered sodium chloride peptone solution pH 7.0. For fatty substances add to the liquid polysorbate 20 or polysorbate 80. Transfer one of the membrane filters, intended for the enumeration of bacteria, to the surface of a plate of casein soybean digest agar and the other, intended for the enumeration of fungi, to the surface of a plate of Sabouraud dextrose agar with antibiotics.

Incubate the plates for 5 days, unless a more reliable count is obtained in shorter time, at 30 °C to 35 °C in the test for bacteria and 20 °C to 25 °C in the test for fungi. Count the number of colonies that are formed. Calculate the number of microorganisms per g or per ml of the extract preparation being examined, if necessary count bacteria and fungi separately.

C-7.6 Plate Count for Bacteria

Using Petri dishes 9 to 10 cm in diameter, add to each dish a mixture of 1 ml of the pretreated extract preparation and about 15 ml of liquified casein soybean digest agar at not more than 45 °C.

Alternatively, spread the pretreated extract preparation on the surface of the solidified medium in a petri dish of the same diameter. If necessary, dilute the pretreated extract preparation as described above so that a colony count of not more than 300 may be expected. Prepare at least two such petri dishes using the same dilution and incubate at 30 °C to 35 °C for 5 days, unless a more reliable count is obtained in a shorter time. Count the number of colonies that are formed. Calculate the results using plates with the greatest number of colonies but taking 300 colonies per plate as the maximum consistent with good evaluation.

C-7.7 Plate Count for Fungi

Proceed as described in the test for bacteria but use Sabouraud dextrose agar with antibiotics in place of casein soybean digest agar and incubate the plates at 20 °C to 25 °C for 5 days, unless a more reliable count is obtained in a shorter time. Calculate the results using plates with not more than 100 colonies.

C-7.8 Multiple Tube or Serial Dilution Method

In each of fourteen test tubes of similar size place 9.0 ml of sterile fluid soybean casein digest medium. Arrange twelve of the tubes in four sets of three tubes each. Put aside one set of three tubes to serve as controls. Into each of three tubes of one set ("100") and into fourth tube (A) pipette 1 ml of the solution of suspension of the test specimen (extract) and mix. From tube A pipette 1 ml of its contents into the one remaining tube (B) not included in the set and mix. These two tubes contain 100 mg (or 100 μl) and 10 mg (or 10 μl) of the specimen respectively. Into each of the second set ("10") of three tubes pipette 1 ml from tube A, and into each tube of the third set ("1") pipette 1 ml from tube B. Discard the unused contents of tube A and B. Close well and incubate all of the tubes.

Following the incubation period, examine the tubes for growth. The three control tubes remain clear. Observations in the tubes containing the test specimen, when interpreted by reference to Table 2, indicate the most probable number of microorganisms per g or per ml of the test specimen.

Table 2 Most probable Total Count by Multiple Tube or Serial Dilution Method
(Clause C-7.8)

SI No	Observed Combination of Numbers of Tubes Showing Growth in Each Set			Most Probable Number of Microorganisms per g or per ml
	Number of mg (or ml) of Specimen per Tube			
	100 (100 µL)	10 (10 µL)	1 (1 µL)	
(1)	(2)	(3)	(4)	(5)
i)	3	3	3	>1100
ii)	3	3	2	1100
iii)	3	3	1	500
iv)	3	3	0	200
v)	3	2	3	290
vi)	3	2	2	210
vii)	3	2	1	150
viii)	3	2	0	90
ix)	3	1	3	160
x)	3	1	2	120
xi)	3	1	1	70
xii)	3	1	0	40
xiii)	3	0	3	95
xiv)	3	0	2	60
xv)	3	0	1	40
xvi)	3	0	0	23

C-8 TESTS FOR SPECIFIED MICROORGANISMS

C-8.1 Pre-treatment of the Extract Sample Being Examined

Proceed as described under the test for total aerobic microbial count but using lactose broth or any other suitable medium shown to have no antimicrobial activity under the conditions of test in place of buffered sodium chloride peptone solution pH 7.0.

C-8.2 *Escherichia coli*

Place the prescribed quantity in a sterile screw capped container, add 50 ml of nutrient broth, shake, allow to stand for 1 h (4 h for gelatin) and shake again. Loosen the cap and incubate at 37 °C for 18 to 24 h.

C-8.2.1 Primary Test

Add 1.0 ml of the enrichment culture to a tube containing 5 ml of MacConkey broth. Incubate in a water bath at 36 °C to 38 °C for 48 h. If the contents of the tube show acid and gas, carry out the secondary test.

C-8.2.2 Secondary Test

Add 0.1 ml of the contents of the tubes containing (a) 5 ml of MacConkey broth and (b) 5 ml of peptone water. Incubate in a water bath at 43.5 °C to 44.5 °C for 24 h and examine tube (a) for acid and gas and tube (b) for indole. To test for indole, add 0.5 ml of Kovac's reagent, shake well and allow to stand for 1 min. If a red colour is produced in the reagent layer indole is present. The presence of acid and gas and of indole in the secondary test indicates the presence of *Escherichia coli*.

WC DRAFT DOCUMENT: FOR COMMENTS ONLY
(Not to be reproduced without permission of BIS or used as Indian Standard)

AYD-04 (26483)
October, 2024

Carry out a control test by repeating the primary and secondary tests, adding 1.0 ml of the enrichment culture and a volume of broth containing 10 to 50 *Escherichia coli* (NCTC 9002) organisms, prepared from a 24 h culture in nutrient broth, to 5 ml of MacConkey broth. The test is not valid unless the results indicate that the control contains *Escherichia coli*.

C-8.2.3 Alternative Test

By means of an inoculating loop, streak a portion from the enrichment culture (obtained in the previous test) on the surface of MacConkey agar medium. Cover and invert the dishes and incubate. Upon examination, if none of the colonies are brick red in colour and have a surrounding zone of precipitated bile the sample meets the requirements of the test for the absence of *Escherichia coli*. If the colonies described above are found, transfer the suspect colonies individually to the surface of Levine eosin methylene blue agar medium, plated on petri dishes. Cover and invert the plates and incubate. Upon examination, if none of the colonies exhibits both a characteristic metallic sheen under reflected light and a blue-black appearance under transmitted light, the sample meets the requirements of the test for the absence of *Escherichia coli*. The presence of *Escherichia coli* may be confirmed by further suitable cultural and biochemical tests.

C-8.3 Salmonella

Transfer a quantity of the pretreated preparation being examined containing 1 g or 1 ml of the product to 100 ml of nutrient broth in a sterile screw capped jar, shake, allow to stand for 4 h and shake again. Loosen the cap and incubate at 35 °C to 37 °C for 24 h.

C-8.3.1 Primary Test

Add 1.0 ml of the enrichment culture to each of the two tubes containing (a) 10 ml of selenite F broth and (b) tetrathionate bile brilliant green broth and incubate at 36 °C to 38 °C for 48 h. From each of these two cultures subculture on at least two of the following four agar media: bismuth sulphate agar, brilliant green agar, deoxycholate citrate agar and xylose lysine deoxycholate agar. Incubate the plates at 36 °C to 38 °C for 18 to 24 h. Upon examination, if none of the colonies conforms to the description given in Table 3, the sample meets the requirements of the test for the absence of the genus *Salmonella*. If any colonies conforming to the description in Table 3 are produced, carry out the secondary test.

Table 3 Interpretation of the Confirmatory Tests for *Salmonella*
(Clause C-8.3.1)

Sl No	Medium	Description of Colony
(1)	(2)	(3)
i)	Bismuth sulphite agar	Black or green
ii)	Brilliant green agar	Small, transparent and colorless, or opaque, pinkish or white (frequently surrounded by a pink or red zone)
iii)	Deoxycholate citrate agar	Colorless and opaque, with or without black centers
iv)	Xylose lysine deoxycholate agar	Red with or without black centres

C-8.3.2 Secondary Test

Subculture any colonies showing the characteristics given in Table 3 in triple sugar iron agar by first inoculating the surface of the slope and then making a stab culture with the same inoculating needle, and at the same time inoculate a tube of urea broth. Incubate at 36 °C to 38 °C for 18 to 24 h. The formation of acid and gas in the stab culture (with or without concomitant blackening) and the absence

WC DRAFT DOCUMENT: FOR COMMENTS ONLY
(Not to be reproduced without permission of BIS or used as Indian Standard)

AYD-04 (26483)
October, 2024

of acidity from the surface growth in the triple sugar iron agar, together with the absence of a red color in the urea broth, indicate the presence of *Salmonella*. If acid but no gas is produced in the stab culture, the identity of the organisms should be confirmed by agglutination tests.

Carry out the control test by repeating the primary and secondary tests using 1.0 ml of the enrichment culture and a volume of broth containing 10 to 50 *Salmonella abony* (NCTC 6017) organisms, prepared from a 24 h culture in nutrient broth, for the inoculation of the tubes (a) and (b). The test is not valid unless the results indicate that the control contains *Salmonella*.

A-8.4 *Pseudomonas aeruginosa*

Pretreat the preparation being examined as described above and inoculate 100 ml of fluid soybean casein digest medium with a quantity of the solution, suspension or emulsion thus obtained containing 1 g or 1 ml of the preparation being examined. Mix and incubate at 35 °C to 37 °C for 24 to 48 h. Examine the medium for growth and if growth is present, streak a portion of the medium on the surface of cetrimide agar medium, each plated on petri dishes. Cover and incubate at 35 °C to 37 °C for 18 to 24 h.

If, upon examination, none of the plates contains colonies having the characteristics listed in Table 4 for the media used, the sample meets the requirement for freedom from *Pseudomonas aeruginosa*. If any colonies conforming to the description in Table 4 are produced, carry out the oxidase and pigment tests. Streak representative suspect colonies from the agar surface of cetrimide agar on the surfaces of *Pseudomonas* agar medium for detection of fluorescein and *Pseudomonas* agar medium for detection of pyocyanin contained in petri dishes. Cover and invert the inoculated media and incubate at 33 °C to 37 °C for not less than 3 days. Examine the streaked surfaces under ultra violet light. Examine the plates to determine whether colonies conforming to the description in Table 4 are present. If growth of suspect colonies occurs, place 2 or 3 drops of a freshly prepared 1 percent w/v solution of N,N,N1,N1-tetramethyl-4-phenylenediamine dihydrochloride on filter paper and smear with the colony. If there is no development of a pink color, changing to purple, the sample meets the requirements of the test for the absence of *Pseudomonas aeruginosa*.

Table 4 Interpretation of the Confirmatory Tests for *Pseudomonas aeruginosa*
(Clause C-8.4)

Sl No	Selective Medium	Characteristic Colonial Morphology	Fluorescence in UV Light	Oxidase Test	Gram Stain
(1)	(2)	(3)	(4)	(5)	(6)
i)	Cetrimide agar	Generally greenish	Greenish	Positive	Negative rods
ii)	<i>Pseudomonas</i> agar medium for detection of fluorescein	Generally colorless to yellowish	Yellowish	Positive	Negative rods
iii)	<i>Pseudomonas</i> agar for detection of pyocyanin	Generally greenish	Blue	Positive	Negative rods

C-8.5 *Staphylococcus aureus*

Proceed as described under *Pseudomonas aeruginosa*, if upon examination of the incubated plates, none of them contains colonies having the characteristics listed in for the media used, the sample meets the requirements for the absence of *Staphylococcus aureus*. If growth occurs, carry out the coagulase test. Transfer representative suspect colonies from the agar surface of any of the media listed in Table 5 to

WC DRAFT DOCUMENT: FOR COMMENTS ONLY
(Not to be reproduced without permission of BIS or used as Indian Standard)

AYD-04 (26483)
October, 2024

individual tubes, each containing 0.5 ml of mammalian, preferably rabbit or horse, plasma with or without additives.

Table 5 Interpretation of the Confirmatory Tests for *Staphylococcus aureus*
(Clause C-8.5)

Sl No.	Selective Medium	Characteristic Colonial Morphology	Gram Stain
(1)	(2)	(3)	(4)
i)	Vogel Johnson agar	Black surrounded by yellow zones	Positive cocci (in clusters)
ii)	Mannitol salt agar	Yellow colonies with yellow zones	Positive cocci (in clusters)
iii)	Baird parker agar	Black, shiny, surrounded by clear zones of 2 to 5 mm	Positive cocci (in clusters)

Incubate in water bath at 37 °C examining the tubes at 3 h and subsequently at suitable intervals up to 24 h. If no coagulation in any degree is observed, the sample meets the requirements of the test for the absence of *Staphylococcus aureus*.

C-8.6 Validity of the Tests for Total Aerobic Microbial Count

Grow the following test strains separately in tubes containing fluid soybean casein digest medium at 30 °C to 35 °C for 18 to 24 h or, for *Candida albicans*, at 20 °C for 48 h.

<i>Staphylococcus aureus</i>	(ATCC 6538; NCTC 10788)
<i>Bacillus subtilis</i>	(ATCC 6633; NCIB 8054)
<i>Escherichia coli</i>	(ATCC 8739; NCIB 8545)
<i>Candida albicans</i>	(ATCC 2091; ATCC 10231)

Dilute portions of each of the cultures using buffered sodium chloride peptone solution pH 7.0 to make test suspensions containing about 100 viable microorganisms per ml. Use the suspension of each of the microorganisms separately as a control of the counting methods, in the presence and absence of the preparation being examined, if necessary.

A count for any of the test organisms differing by not more than a factor of 10 from the calculated value for the inoculum should be obtained. To test the sterility of the medium and of the diluent and the aseptic performance of the test, carry out the total aerobic microbial count method using sterile buffered sodium chloride peptone solution pH 7.0 as the test preparation. There should be no growth of microorganisms.

C-8.7 Validity of the Tests for Specified Microorganisms

Grow separately the test strains of *Staphylococcus aureus* and *Pseudomonas aeruginosa* in fluid soybean casein digest medium and *Escherichia coli* and *Salmonella typhimurium* at 30 °C to 35 °C for 18 to 24 h. Dilute portions of each of the cultures using buffered sodium chloride peptone solution pH 7.0 to make test suspensions containing about 10³ viable microorganisms per ml. Mix equal volume of each suspension and use 0.4 ml (approximately 10² microorganisms of each strain) as an inoculum in the test for *E. coli*, *S. typhimurium*, *P. aeruginosa* and *S. aureus*, in the presence and absence of the extract preparation being examined, if necessary. A positive result for the respective strain of microorganism should be obtained.